



commode assembly

1. Remove contents from box. Place side frame on the floor with the receivers facing up. Install the center frame into the side frame by depressing the brass pin and inserting both the front and rear pins into the side frame. The center frame should be inserted first into the receiver that runs parallel to the side frame (this is also the front side of your commode). Once this is locked into position, take the center frame and insert it into the receiver that runs perpendicular to the side frame. (This is the back side of your commode). Make sure the brass pins clearly snap into position by pulling back to insure the pins are fully locked into position.

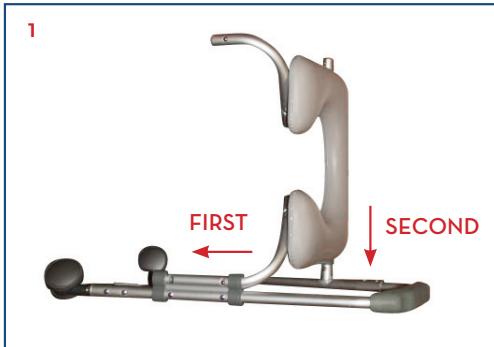


2. Place the second side frame on the floor with the receivers facing up. Flip the side frame which is attached to the center section and install the other side of the center section into the side frame by depressing the brass pin and inserting both the front and rear pins into the side frame. The same procedure should be followed by locking the front of the commode in first and then locking the back of the commode into place. Make sure the brass pins clearly snap into position by pulling back to insure the pins are fully locked into position.

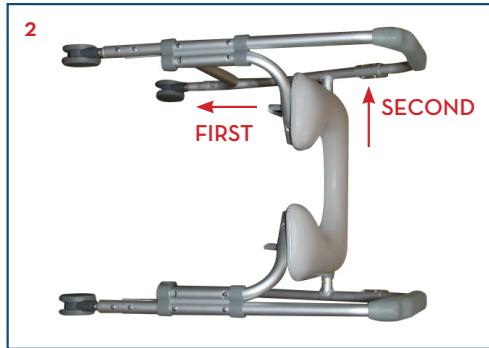

3. To install the backrest, insert backrest into backrest receptacles by depressing the brass pins on both sides into the commodes back frame. The brass push pins will snap into place when secure. Check by pulling back to insure the back is secure and locked into place.



4. To install the front foot bar depress the the brass button and insert the front foot bar and install into the front receiving holes. You will need to depress the brass button to insert the bar into the receiving tube and make sure the brass pins clearly snap into place.

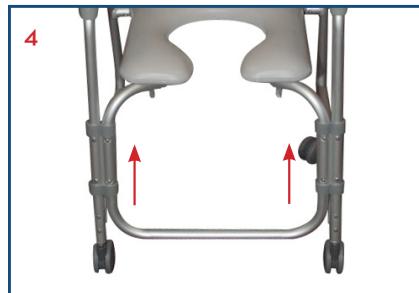

5. Your commode comes with a bucket and lid as well as a splash shield. The bucket is used while the commode is used by the bedside. The splash shield is used while the commode is used over the toilet. To install position the bucket so two of the bucket tabs are on the bucket rails on the back side of the commode. Gently slide the bucket through the rails until the bucket is properly positioned on the front of the commode. To install the splash shield place the commode over the toilet and place the splash shield so two of the splash shield tabs are on the bucket rail.



additional precautions

1. Do NOT install commode without first reading and understanding this assembly and instruction sheet. If you can not understand these instructions, contact a healthcare, dealer or technical professional before attempting to install. Otherwise, injury or damage may occur.
2. Inspect all parts for shipping damage before assembly. If there is shipping damage - DO NOT USE. Contact dealer/carrier for further instructions.
3. Ensure ALL screws, nuts and/or bolts are tightened.
4. All four legs MUST touch the floor simultaneously at all times.
5. Users with limited physical strength should be supervised or assisted while using commode.