

3M™ Tegaderm™ I.V. Advanced Securement Dressings
Product Information



Designed for Securement & Comfort

3M™ Tegaderm™ I.V. Advanced Securement Dressings are comprehensively designed to deliver exceptional patient care through advanced catheter securement, superior wear time, gentle removal, and excellent moisture management.

- A line of dressings which optimizes all characteristics to provide advanced catheter securement
- A cost-effective securement solution
- Innovative adhesives balance securement and gentle removal
- Waterproof and highly breathable to provide better wear time
- Easy to apply and remove

Tegaderm™ I.V. Advanced
Part of the Full Line of Tegaderm™ Brand I.V. Site Dressings



3M™ Tegaderm™ I.V. Advanced Securement Dressings

Provides secure hold

- Reinforced border, deep notch, and tape strips work together to enlist the entire dressing in securement
- Innovative adhesives hold strongly, manage moisture, and release gently
- Designed to maintain securement by preventing edge lift, flexing with patient movement, and managing moisture

Excellent value


- Can potentially reduce the number of dressing changes and restarts
- May be worn for the life of a PIV catheter
- Provides up to 7 days of wear time for CVC's**

Easy to apply and remove

- Designed to promote consistent application
- Frame delivery makes placement accurate and easy

Supports infection prevention goals

- A waterproof film coating on all pieces provides a barrier to external contaminants including liquids, bacteria and viruses*
- Allows continuous monitoring

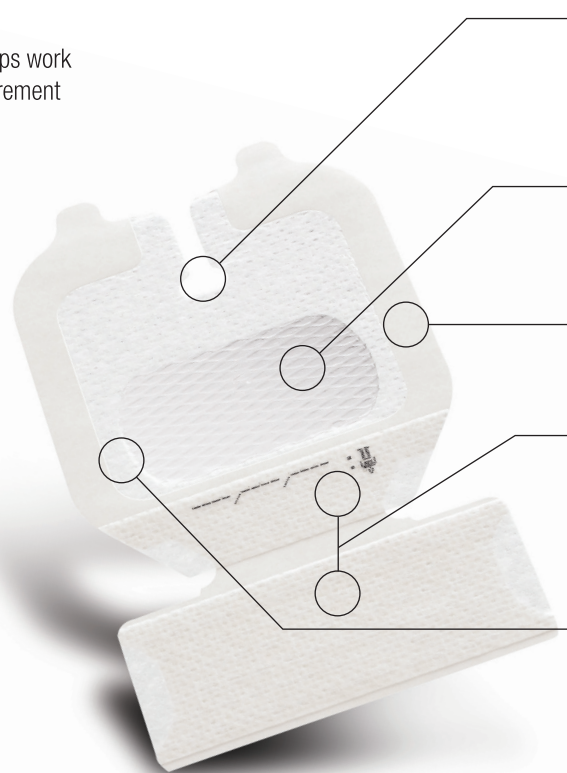

Deep notch

- Fits easily and conforms around hubs and lumens
- Designed to discourage edge lift
- Designed to reduce the risk of catheter manipulation and dislodgement

Highly breathable transparent film

- Designed to release moisture quickly to provide better wear time and securement
- A waterproof, sterile barrier to external contaminants*

Picture-frame delivery

- Minimizes potential to stick to gloves or to itself
- One-handed application

Securement tape strips

- Enhance securement while promoting consistent application
- Waterproof film coating
- Pre-printed documentation label for dressing changes

Stabilization border

- Maximizes securement, breathability and wear time
- Flexible and conformable
- Reduces edge lift to maintain securement
- Waterproof film coating

* *In vitro* testing shows that the transparent film of 3M™ Tegaderm™ I.V. Advanced Securement Dressing provides a viral barrier for viruses 27 nm in diameter or larger while the dressing remains intact without leakage.

** O'Grady, N et al. (2011). Guidelines for the Prevention of Intravascular Catheter-Related Infections, 2011. Retrieved from <http://www.cdc.gov/hicpac/BSI/BSFguidelines-2011.htm>. (pp. 32-35)

Product and package are latex free.

Ordering Information

PRODUCT CODE	OVERALL DRESSING SIZE	SUGGESTED DEVICES	DRESSINGS/ BOX	BOXES/ CASE
1683	2½ in. x 2¾ in. (6.5 cm x 7 cm)	PIV	100	4
1685	3½ in. x 4½ in. (8.5 cm x 11.5 cm)	CVC, PICC	50	4

3M™ Tegaderm™ Dressings SIMPLE. DEPENDABLE. TRUSTED.



Medical Solutions Division

3M India Limited

24, Concorde Block, UB City,
Vittal Mallya Road Bangalore I 560001
E-mail: 3mhealthcare.india@mmm.com
Toll Free 1-800-425-3030

Printed on MARCH-2020