

Delta-Dry®

Water Resistant Cast Padding



The solution that fits your life

- Conformable
- Water resistant
- Traditional cast removal with no extra components needed



Delta-Dry®

Water Resistant Cast Padding

Delta-Dry®'s unique patented technology outperforms other water resistant cast padding*, making Delta-Dry® the best choice for medical professionals and a more convenient solution for their patients over traditional cast padding.

* refer to chart below

Product Benefits

- Allows patients to shower, swim and bathe like normal
- No special tools are needed for application or removal
- Conformable enabling patient to having increased comfort and protection
- Easy to apply
- Better long term water resistance
- Offers more under cast protection
- Eliminates the need for stockinette




Directions

- The open knit structure, along with the water resistant qualities of the fibers, allows water to quickly drain from the cast.
 - Position the cast to allow water to drain from the cast for approximately 10 minutes.
 - Long arm and leg casts must be positioned as to not let water remain in the heel or elbow.

Specifications for Delta-Dry®

Catalogue Number	Dimensions	Box Contents
73443-00	5cm x 2.4m (2" x 2.6 yds)	12
73443-01	7.5cm x 2.4m (3" x 2.6 yds)	12
73443-02	10cm x 2.4m (4" x 2.6 yds)	12

MVTR Delta-Dry® vs. Brand X

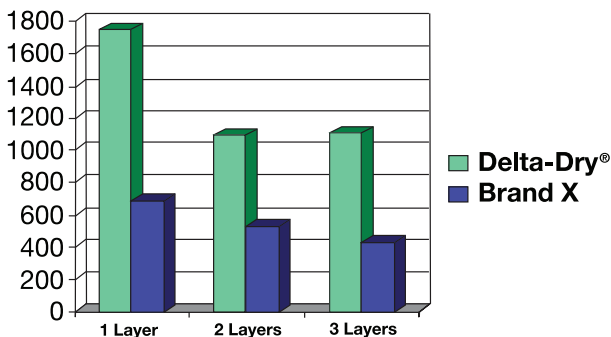

Moisture Vapour Transmission Rate

- Independent laboratory testing compared the Moisture Vapour Transmission Rate (MVTR)

Laboratory testing results are intended to illustrate a product's performance under controlled conditions. Actual use results may vary.

