

WALKER AMBULATORY SUPPORT



«The first step is the hardest.» This old proverb suddenly seems more relevant than ever for residents facing a new challenge: to regain their ability to walk.

The Walker® provides the optimum conditions to take that first step – a safe environment for both resident and therapist. And the Walker goes on to play an important supporting role every step of the way during the initial stages of a resident's rehabilitation and ambulation training.

The Walker's hi-low function reduces transfer risks for both resident and physiotherapist when hoisting or lowering the resident to and from the standing position. While in motion, the resident's sense of security is further enhanced by safety and comfort features: adjustable handgrips for a firm hold, the U-shaped support table for leaning on and back/suspension straps to prevent falling.



...with people in mind

FLEXIBILITY FOR THE DIVERSE CHALLENGES OF EVERYDAY CARE



Once the resident is in a comfortable standing position, the spade seat can be easily removed. (For the sake of clarity, here shown without resident).



Adjustable handgrips offer the resident a firm hold.




A straight steering device makes the Walker even sturdier, enabling the resident to walk straight forward. The Walker is equipped with brakes on two castors.



The battery-powered version of Walker is equipped with an easy-to-use hand control that allows the therapist to stay close to her resident.



There is no need to park the Walker to recharge it – you simply exchange the battery for a fully charged one from the battery charger.

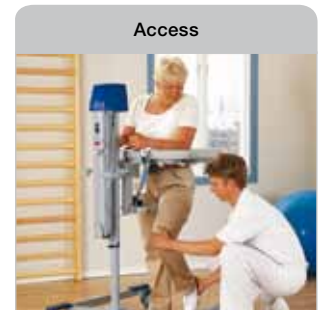

The Walker's open sides offer the therapist free access to the resident's hips, knees and ankles.



Suspension straps, together with the back strap, mean there is no risk that the resident can fall.



The U-shaped support table is soft and comfortable to lean on.

Product information

Safe working load 136 kg (300 lbs)

Lifting stroke range (battery powered) 565 mm (23 1/4")
(hydraulic) 500 mm (19 3/4")

Width 825 mm (32 1/2")

Max length (hydraulic) 1040 mm (41")
Max length (electric) 1000 mm (39 1/2")

Battery (removable & rechargeable) 24V DC

Table height (hydraulic) 770-1270 mm (30 1/4 - 50")
(electric) 810-1375 mm (32 - 54 1/4")

Spade seat height (hydraulic) 360-860 mm (14 1/4 - 33 3/4")
(electric) 400-965 mm (15 3/4 - 38")

Weight (hydraulic) 42.5 kg (93.5 lbs)
(electric) 48.5 kg (107 lbs)

Only ArjoHuntleigh designed parts, which are designed specifically for the purpose, should be used on the equipment and products supplied by ArjoHuntleigh. As our policy is one of continuous development we reserve the right to modify designs and specifications without prior notice.

© and ™ are trademarks belonging to the ArjoHuntleigh group of companies. © ArjoHuntleigh, 2011.