




DONJOY®




ICEMAN® CLEAR³

Individual results may vary. Neither DJO Global, Inc. nor any of its subsidiaries dispense medical advice. The contents of this brochure do not constitute medical, legal, or any other type of professional advice. Rather, please consult your healthcare professional for information on the courses of treatment, if any, which may be appropriate for you.

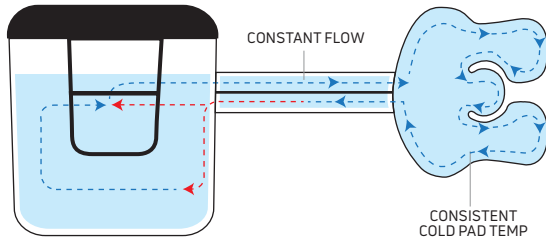

ICEMAN with dual pump recirculation technology. The DonJoy® IceMan® features a recirculation system that allows water warmed after flowing through the cold pad to be preserved and remixed with cooler ice water at a constant flow rate, providing consistent cool water distribution throughout the cold pad.

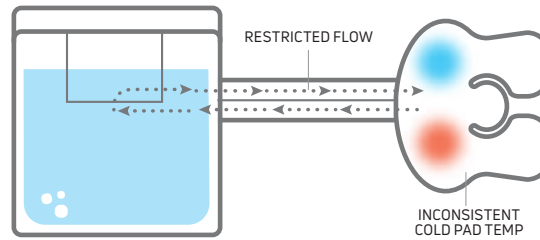

COMPETITORS Other versions restrict or pinch water flow to the cold pad, causing variable flow rates and inconsistent cold pad temperatures.



WARNING!

THIS DEVICE CAN BE COLD ENOUGH TO CAUSE SERIOUS INJURY. Please read and understand all warnings and instructions for use before using the device. DO NOT use this device without a prescription and direction from a physician.

DONJOY® ICEMAN® CLEAR³

The IceMan® CLEAR³ cold therapy unit helps reduce pain and swelling, speeding up rehabilitation and recovery. The IceMan helps provide extended cold therapy to aid in a variety of indications and protocols as directed by a medical professional. It utilizes DonJoy's recirculation system, which helps maintain more consistent and accurate temperatures than other competitive cold therapy units, in a pre-set configuration.

PROVEN PERFORMANCE

- Time tested with over 20-years of experience delivering industry leading cold therapy solutions
- Developed through direct customer interaction and feedback

TECHNOLOGY ADVANTAGE

- Recirculation system helps deliver consistent cold therapy throughout the cold pad, eliminating freezing inlet temperatures

PATIENT CONVENIENCE AND EASE-OF-USE

- Clear cooler design for easy refill indication
- Self-priming, plug-n-play operation
- Designed for easy portability and storage
- A complete line of cold pads to help deliver cold therapy to where it is needed the most

The use of cold therapy may reduce the need for narcotics and help accelerate rehabilitation and recovery.

STEPS:

1. Add ice to fill line inside the device.
2. Add cold water to fill line.
3. Place lid on the device making sure the lip inserts in groove. Then press front of the lid down to close and secure.
4. Connect the IceMan® CLEAR³ hose to the cold pad hose. To ensure a reliable connection, "snap" or "click" hoses together into place so that the fit is tight and snug.
5. To turn device on, insert cord into the connection port on the back of the device and plug the power supply into the wall outlet.