


Seating - Corner Chairs

Kaye Corner Chairs

Corner Chairs are for users who need support of the head, trunk, and pelvis in order to both develop an upright sitting posture and align the head and trunk. Two sizes are available, large and small.

Kaye Corner Chairs have a firm, padded insert that converts the V-shaped back into a three-sided winged back, providing support along the spine, giving the child better spinal extension and resulting in greater freedom of the head and eyes for scanning and attending. This insert can be removed so that the chair sides provide even greater assistance in bringing the arms forward on the tray.

The seats on the Kaye Corner Chairs can be raised so a child can be positioned with the hips and knees flexed or it can be used in a lower, long-sitting position.

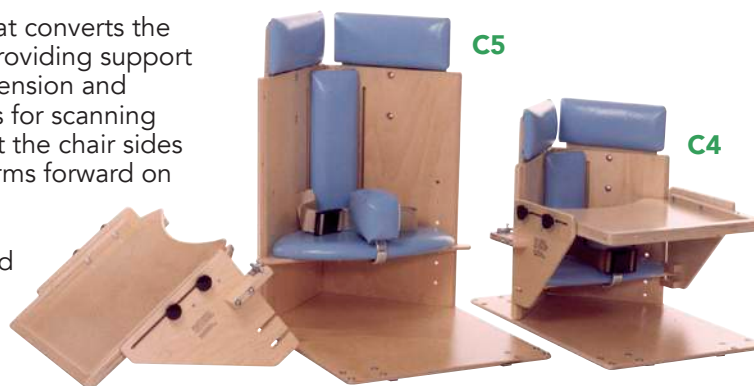

Model Features	C4	C5
Age	8 months to 4 years	4 to 8 years
Size & Weight	Up to 39" / 40 lbs (18 kg)	Up to 58" / 65 lbs (29 kg)
Height of chair with head supports	Adjusts from 20" to 27"	Adjusts from 29" to 41"
Removeable Tray	15" x 17"	18" x 20"
Seat depth with insert	8"	10.5"
Seat depth without insert	12"	14"
Removable abductor	3" wide	3" wide
Chair with 4" swivel locking casters	Model C4C	Model C5C

