



# F5 Corpus



Fully independent oil-dampened suspension for a more stable, comfortable ride



Lower 17.5" seat-to-floor height for improved accessibility

**THE F5 BASE.** It features an array of advancements to help increase the stability of the chair for a smooth, comfortable ride, including: a smaller footprint with an improved center of gravity, a more rigid chassis, fully independent oil-dampened suspension and a lower seat-to-floor height. The new three-point tilt and elevator system provides solid stability and performance. The standard 7.5 mph motor package with improved torque allows users to easily navigate obstacles at any speed. Options include: black tires with more aggressive tread patterns, chrome wheel accents and two additional color options.

**THE CORPUS® SEATING SYSTEM.** We retained all of the things that you have come to expect from the legendary Corpus seating system like the arc legrests, the open armrest design and the Uni-Track System. Then, we took the Corpus seating system to the next level with advanced features such as: the 14" seat elevator and the A/P tilt function, which allows for optional 5°, 10°, 20° and 45° of anterior tilt, while maintaining 50° of posterior tilt. The A/P tilt helps our users with transfers, added reach, and allows them to get closer to a table or workspace.



Enhanced tilt functionality  
with optional 5°, 10°, 20° and 45° anterior tilt, while maintaining 50° posterior tilt

Anterior tilt must be clinically recommended by a licensed professional. Additional support features such as knee blocks and chest straps may need to be prescribed.



Power everything  
including tilt, recline, legrests and more

14" seat elevator  
that maintains three-points of  
contact for increased stability



**permobil**

THE POWER OF MOBILITY™




**F5**  
**Corpus**

FEATURING A/P TILT

# F5 Corpus



## AVAILABLE FABRIC

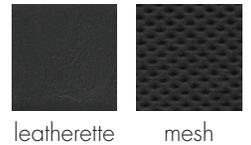

leatherette

mesh

## AVAILABLE TRIM COLORS



galactic green



midnight black



ocean blue



phantom grey



power pink



radiant red



sunburst orange



titan blue

## SPECIFICATIONS

Maximum User Weight	300 lbs
Maximum Speed	7.5 mph 6.0 mph without anti-tippers
Range	16 mi*
Base Width	25.5"
Base Length	42" with anti-tippers 38" without anti-tippers
Minimum Turning Radius	30"
Minimum Ground Clearance	3"
Recommended Battery Type	Group 24
Weight, including batteries	410 lbs
Drive Electronics	R-Net 120A

Seat-to-Floor Height (STFH)	17.5", 18.5", 19.5"
Optional Seat Elevator	14"
Tilt Options	0° - 50° Posterior 5°, 10°, 20° & 45° Anterior
Recline Options	85° - 120° (manual) 85° - 180° (power)
Legrest Angle	85° - 170°
Backrest Height	20", 23" - 28" (by 1" incr.)
Seat Depth	14" - 22" (by 1" incr.)
Seat Width	17", 19", 21", 23"
Armrest Pad Lengths	10", 13", 16", 18"

\* Actual driving range will vary based on driving and battery conditions.