



caribou

by R82[®]



New Caribou standing frame

R82 A/S

Tel. +45 79 68 58 88 - R82@R82.com - www.R82.com

R82 UK Ltd.

Tel. 0121 561 2222 - R82uk@R82.com - www.R82-uk.co.uk

 nothing compares to a smile!



caribou



R82

caribou[™]

Vertical-horizontal tilt

The Caribou tilts from vertical to horizontal. The horizontal position makes it easier for the carer and less stressful for the user to be placed into position.



Adjustable plate

Height-adjustable for correct positioning.

Shaping of the base plate

Side slot wedges for shaping of the base plate.



Parted foot plates

The foot plates allow positioning and support from almost any angle.



Castors with locks


Secure the Caribou for safe transfers.



The Caribou is a new standing frame. The frame supports the prone or supine position. The Caribou can move from a horizontal to a vertical position and is designed to assist standing in people with a wide range of conditions, including spasticity, contractures and circulatory problems.

The Caribou is built with a base board. A series of optional top supports are available to provide the optimum freedom of movement. A wide range of accessories are also available to meet the specific needs of the user.

The Caribou is available in four sizes and suitable for children between the ages of 1-18.

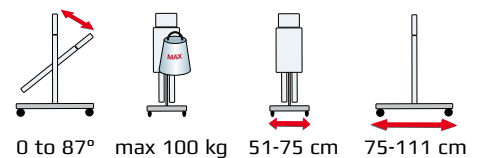



caribou

Accessories

The Caribou comes with a variety of accessories to support the individual child and his or her specific needs.



Side supports



Strap for side supports



Shoulder support



Wedges for side slots



Head plate



Extra cushion



Back support



Pommel



Strap for pommel



Knee support, fixed



Knee support, swing away



Formed foot plates, set



Heel support for parted foot support, set



Tray with bowl and cover





Arm rests for tray



Toybar for tray



Grip handle for tray

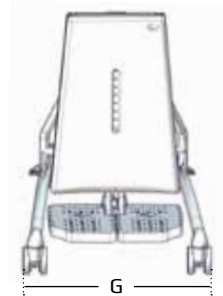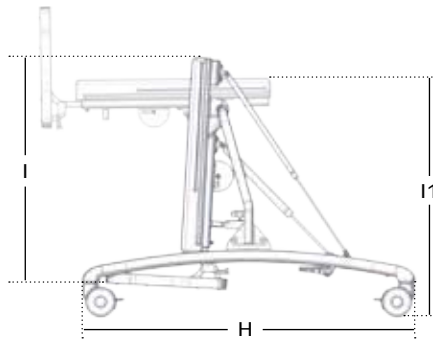

Measurements

	Size 1 cm (inch)	Size 2 cm (inch)	Size 3 cm (inch)	Size 4 cm (inch)
Frame width (G)	51 (20")	57 (22½")	65 (25¼")	75 (29¼")
Frame length (H)	75 (29¼")	92 (36")	111 (43¼")	111 (43¼")
Height, vertical (I)	46-57 (18-22½")	53-64 (20¾-25¼")	68-79 (26¾-31")	88-99 (34¾-39")
Height, horizontal (I1)	70 (27½")	70 (27½")	80 (31½")	80 (31½")
Height from floor to frame	17 (6¾")	18 (7")	18 (7")	18 (7")
Basis plate width, top-bottom	25-32 (9¾-12½")	30-38 (11¾-15")	35-46 (13¾-18")	40-55 (15¾-21¾")
Basis plate height	43 (17")	50 (19¾")	65 (25½")	85 (33½")
Angling	0-87°	0-87°	0-87°	0-87°
Foot plate angling	-10 to 10°	-10 to 10°	-10 to 10°	-10 to 10°
Weight wo/ power	20.5 kg (45.2lb)	24 kg (53lb)	30.5 kg (67.2lb)	35 kg (77.2lb)
Weight w/ power	22.5 kg (49.6lb)	25.5kg (56.2lb)	32.5 kg (71.7lb)	37 kg (81.6 lb)
Max. load	40 kg (88.2lb)	60 kg (132.2lb)	80 kg (176.4lb)	100 kg (220.5lb)