



21-42 cm



18-50 cm



28-50 cm



max. 70 kg

wombat living



Wombat Living exudes Danish design and the minimalist appearance make it fit into any home. It is available in 3 sizes and configurable, making it fit any shape and size.



Home environment • School • Institution

R82[®]

nothing compares to a smile!

+ Danish design | Innovation | Quality

CONTENTS

- 4 When design meets function
- 6 Numerous configurations
- 8 Design your chair
- 10 It is not just about sitting ...
- 11 Covers and colours
- 12 Girls just have more fun
- 14 Lean forward and be active
- 16 Achieve new heights
- 18 Features and functions
- 20 Accessories
- 22 How to measure the user
- 23 Measurements

The choice is yours...

At R82 we spend considerable resources precisely on the development process with the clear purpose of: creating user-friendly products providing pleasure and reserves of energy, and

greater physical well-being. Also we are not satisfied until there is complete harmony between function and design. We recommend Wombat Living as your next choice!

WOMBAT LIVING IS THE ACTIVITY
CHAIR TO USE ALL DAY, EVERYDAY
– AND FOR A LIFE TIME...

R82 





When design meets function!

Through continued development and customer and therapist feedback, the new Wombat Living provides all features and benefits needed for a functional activity chair.


Wombat Living is a cool activity chair for both children and teens. It has built-in growth making it a chair for life. The new fresh look with modern colours makes it suitable for use in any indoor environment.

The height adjustable frame and the anterior tilt function ease standing transfer and ensure that all unnecessary and repetitive lifting is avoided.


Wombat Living exudes Danish design and the minimalist appearance with the discreet frame and unique fabric colours, make it fit anywhere. It is available in 3 sizes and configurable, making it fit any shape and size.



DANISH DESIGN • INNOVATION • QUALITY

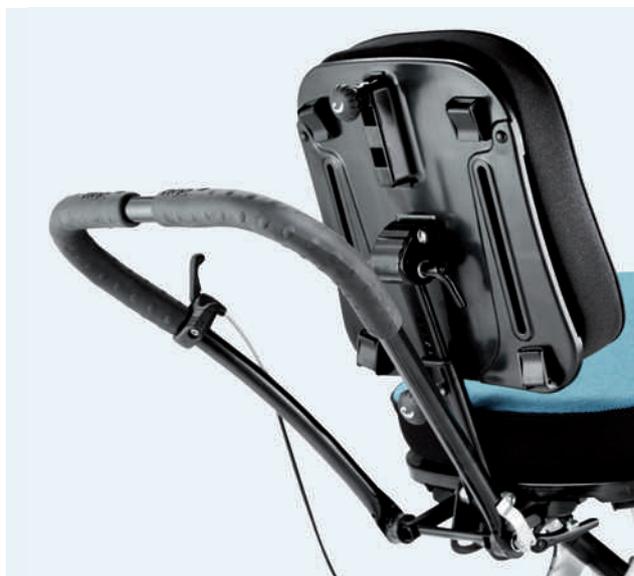

The thoughtful design and mounting of it integrate the new angle adjustable push brace with the chair. Use the push brace to assist the height and angle adjustment

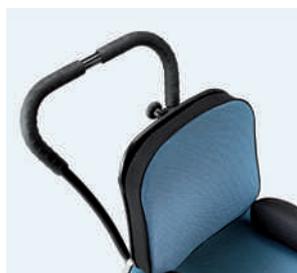

and angle it 'up and away' when it is not in use. Push handles and push brace for mounting on the back of the chair are available as accessories.



Wombat Living is equipped with either four 75 mm or 100 mm smooth running castors, making the chair very easy to manoeuvre. All castors are equipped with



brakes to secure a stable transfer in and out of the chair. A central brake system can be added for easy and quick handling.



A small selection of possibilities...

Numerous configurations

Wombat Living is an activity chair with numerous configuration possibilities. It is perfect for daily activities such as reading and learning at school, having fun with friends and when helping with daily chores at home.

Depending on the user's level of function, the chair can be designed with standard functions that

suit the specific needs. For some users a chair with only a few accessories and simple standard features is suitable. Others might need a chair with much more accessories and technical functions, like a power unit for height adjustment.

Nothing compares to a custom designed chair!



Mix and match!

- 1 Choose size
- 2 Choose frame colour
- 3 Choose castor size
- 4 Choose how to height adjust
- 5 Choose how to angle adjust
- 6 Choose cushion colour and shape



Go to next page



Design your chair...

It is very simple to design your own custom chair. Choose size, frame colour, castors, how to perform the height and seat angle adjustments, and choose your favourite among the selection

of comfortable cushions in new unique fabric colours. No need for further customization or special adjustments. Wombat Living is always delivered - ready for use!



Size 1

Size 2

Size 3



White

Anthracite grey



75 mm



100 mm

Central brake for 100 mm castors

4

HEIGHT ADJUSTMENT



Pedal mounted on the frame



Release handle mounted on the push brace



Power actuator

5

SEAT ANGLE ADJUSTMENT



Release handle mounted beneath the seat



Release handle mounted on the push brace

6

CUSHIONS & COLOURS



Red standard cushions



Petroleum standard cushions



Grey standard cushions



Grey contoured cushions



It is not just about sitting...

If you appreciate good seating comfort as much as we do - Wombat Living is the chair for you!

To meet various postural needs, the cushions are available in different sizes and shapes. The standard cushions provide good seating comfort, whilst the contoured cushions combine comfort

with excellent support for both seat and back. For added support a lumbar support is available as an accessory. The cushions are made from flame-retardant foam meeting all standard requirements.

When seated – it is all about sitting comfortably. Therefore, Wombat Living is your best next choice.



A good seating position must:

- promote correct posture
- prevent pressure points
- be comfortable
- be supportive

Covers and colours

The cover is available in three appealing colours and is breathable and water-repellent for added comfort. Both the seat and back covers are treated to repel dry soil and unpredictable spills and splashes. Any liquid spills can be wiped off and dry

soil can be brushed off easily. Now you can live, work or play with less worry.

A washable incontinence seat cover is available as an accessory.

Small details like cushion shapes that protect the back of the user's knees and extra durable edges make a big difference in everyday use





Lean forward and be active

The Wombat Living has up to -15° forward tilt, a huge benefit for transfer of users with lower limb innervation and whose feet are supported on the ground or assisted for independent transfer.

For more users the active position helps maintain appropriate posture, making it possible to interact even more. Further, studies have shown that the forward sloping chair provides better alignment for the spine and reduces the muscular load on the back.



Achieve new heights

At configuration you choose between three different ways to perform the height adjustment. The pedal mounted on the frame and the release handle mounted on the push brace can be activated by the assistant, while the power

actuator can be activated both by the assistant and the user. The Wombat Living adjusts from seat height 36 cm to 70 cm above floor. The height adjustable frame eases standing transfer into the chair and helps ensure that unnecessary and repetitive lifting is avoided.





FEATURES & FUNCTIONS

1. Comfortable cushions with extra durable edges
2. Unique fabric colour selection
3. Seat depth adjustment for built-in growth
4. Quick and easy tilt-in-space
5. Fix lock options to fit the unique vests & belt programme
6. Prepared for mounting accessories
7. Angle adjustable push brace
8. Height adjustable frame for easy transfer
9. Foot support with adjustable height, depth and angle
10. Smooth-running lockable castors





Accessories

If you're looking to add comfort and functionality to your Wombat Living, a comprehensive range of accessories to meet your needs, is available.

Adding the needed support to the upper trunk, head and arms, the chair lets you rest comfortably.

All accessories always come complete with installation hardware making setup easy.

Upgrade your Wombat Living by adding all sorts of useful functionality with the right accessories.



Checklist!

- Head support
- Back wedge
- Side supports
- Incontinence seat cover
- Lower back extension
- Helper's handle
- Adjustable push brace
- Hip supports
- Pommel
- Arm rests
- Tray
- Anti-tip device




[01] Head support. [02] Back wedge. [03] Side support, fixed. [04] Side supports, swing-away. [05] Anti-tip device. [06] Incontinence seat cover. [07] Lower back extension. [08] Helper's handles. [09] Adjustable push brace. [10] Tray, smoke-coloured. [11] Hip supports. [12] Arm rests. [13] Pommel. [14] Knee abduction pommel.




How to measure the user


This is a simple guide to make sure the chair is correctly configured for the user. To measure accurately, always make sure that the user is lying on a firm flat surface with hips, knees and ankles at a 90 degree angle.

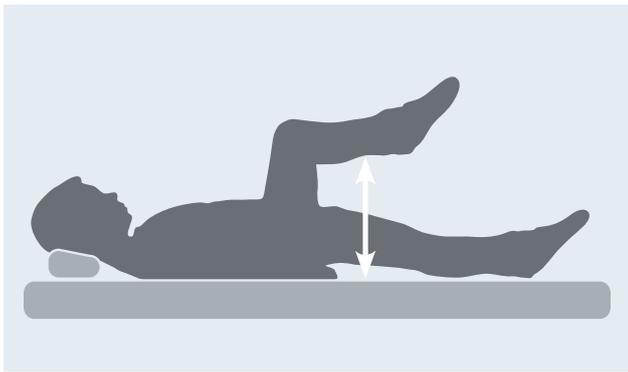


Seat depth

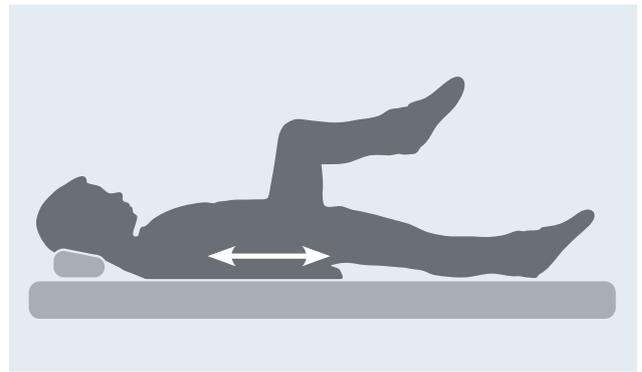


Armpit height

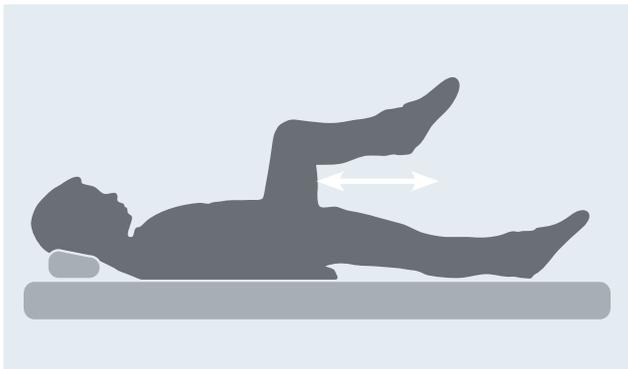

Lower leg length

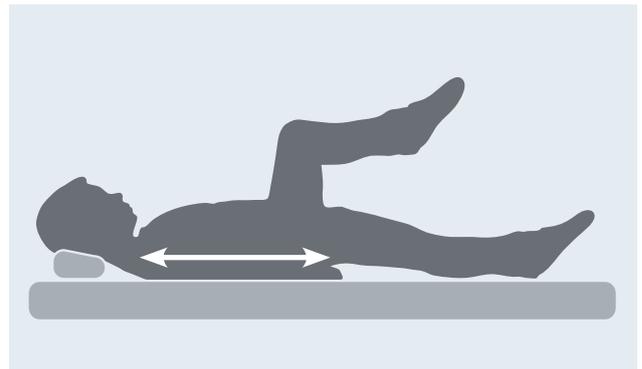


Back height

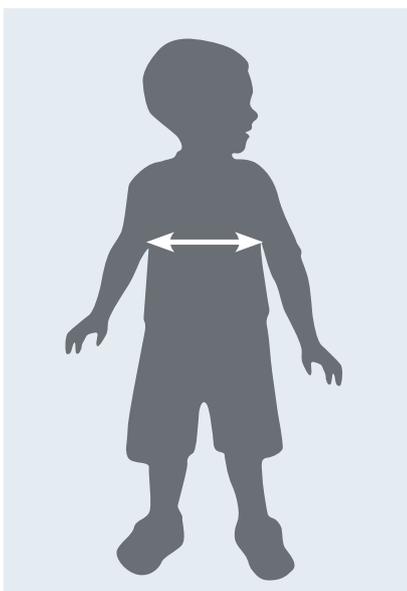

Chest width




Wombat Living is designed to adjust for growth and changing needs. This picture shows how the seat can be altered and illustrates how important good measurement is.

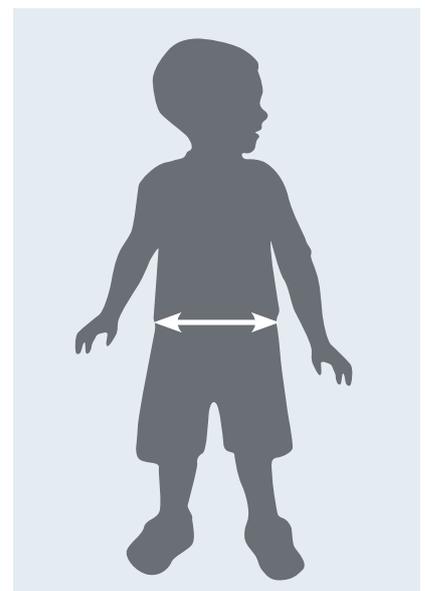

Hip width

Measurements

Wombat Living has built-in growth, meaning adjustable seat depth and hip supports to adjust seat width. Once you have measured the user, you can easily identify the correct chair size in below table.

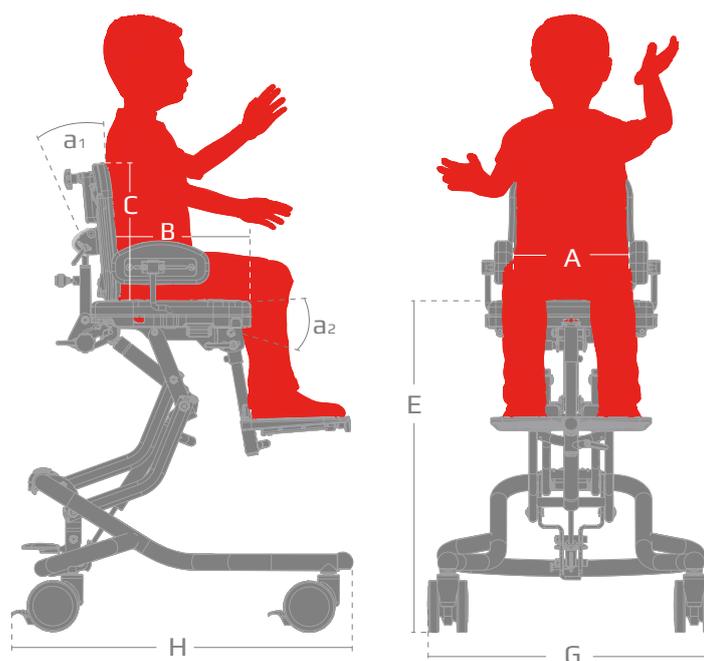

Measurements

	Wombat Living size 1	Wombat Living size 2	Wombat Living size 3
Seat width (A)	30 cm (11¾")	37 cm (14½")	42 cm (16½")
Seat width between hip supports (accessories)	21-30 cm (8¼-11¾")	28-37 cm (11-14½")	25-42 cm (9¾-17¼")
Seat depth (B)	18-31 cm (7-12¼")	25-38 cm (9¾-15")	32-50 cm (12½-19¾")
Back height (C)	28-34 cm (11-13½")	31-38 cm (12¼-15")	39-50 cm (15¼-19¾")
Lower leg length (D)	15-25 cm (6-10")	15-35 cm (6-13¾")	15-44 cm (6-17¼")
Seat height above floor, 75 mm castors (E)	36-65 cm (14¼-25½")	36-65 cm (14¼-25½")	39-67 cm (15¼-26½")
Seat height above floor, 100 mm castors (E)	39-68 cm (15¼-26¾")	39-68 cm (15¼-26½")	42-70 cm (16½-27½")
Back width, flat cushion (F)	25 cm (9¾")	30 cm (11¾")	35 cm (13")
Back width, contoured cushion, between side supports (F)	17 cm (6¾")	23 cm (9")	27 cm (10½")
Arm rest height above seat cushion	6-17 cm (2¼-6¾")	6-17 cm (2¼-6¾")	6-17 cm (2¼-6¾")
Back angle (a1)	90° to 130°	90° to 130°	90° to 133°
Seat angle (a2)	-15° to 20°	-15° to 20°	-10° to 22°
Angle, foot support bar (a3)	-30° to 30°	-30° to 30°	-30° to 30°
Frame width (G)	55 cm (21½")	55 cm (21½")	61 cm (23¾")
Frame length (H)	60 cm (23½")	60 cm (23½")	72 cm (28")
Max. user weight / load	45 kg (99 lb)	45 kg (99 lb)	70 kg (154 lb)
Weight	15,5 kg (34 lb)	15,5 kg (34 lb)	22 kg (48.5 lb)