

600 Mattress

Maximizing support as well as comfort

Pressure support surface with gel-infused open cell memory foam in head and torso sections to provide ample support while keeping the resident comfortable

Advanced pressure support, delivering comfort and relief

Superior pressure redistribution, delivered by enhanced support - means less cradling at higher levels of compression and eliminating static fatigue and shear. Protects vulnerable heel area with 30-degree slope - redistributing pressure to the entire leg, away from the heel. Firm side bolsters provide additional support for easy patient transfer.

Provide great clinical support and comfort for your residents



Measurements:

- Dimensions: 36" x 80"
- Weight capacity: 450 lbs

Materials & Certification:

- High-performing foam: Yes
- Cover: 4 way stretch nylon top cover reduces perspiration and provides a low shear friction surface, that is anti-microbial and fluid resistant

- Bottom cover: Triple laminate vinyl, fluid resistant, anti-microbial/anti-bacterial, fire retardant

- Cover color: Blue
- Vapor-permeable cover: Yes
- Zipper: 360-degree zipper for easy removal, cleaning or replacement. Anti-Contamination Flap covers the zipper, reducing the possibility of fluid migration

- Waterproof: Yes
- Fire retardancy: Fire sleeve meets Federal Fire Code 16 CFR part 1633

Warranty:

5 years for foam, 2 years for cover

1. Side Bolsters

Firm high density foam

2. Base layer

Firm high density foam circulation cuts to prevent tearing

3. Middle layer


Firm medium density foam

4. Top layer

Gel infused memory foam

5. Heel Slope

Soft foam

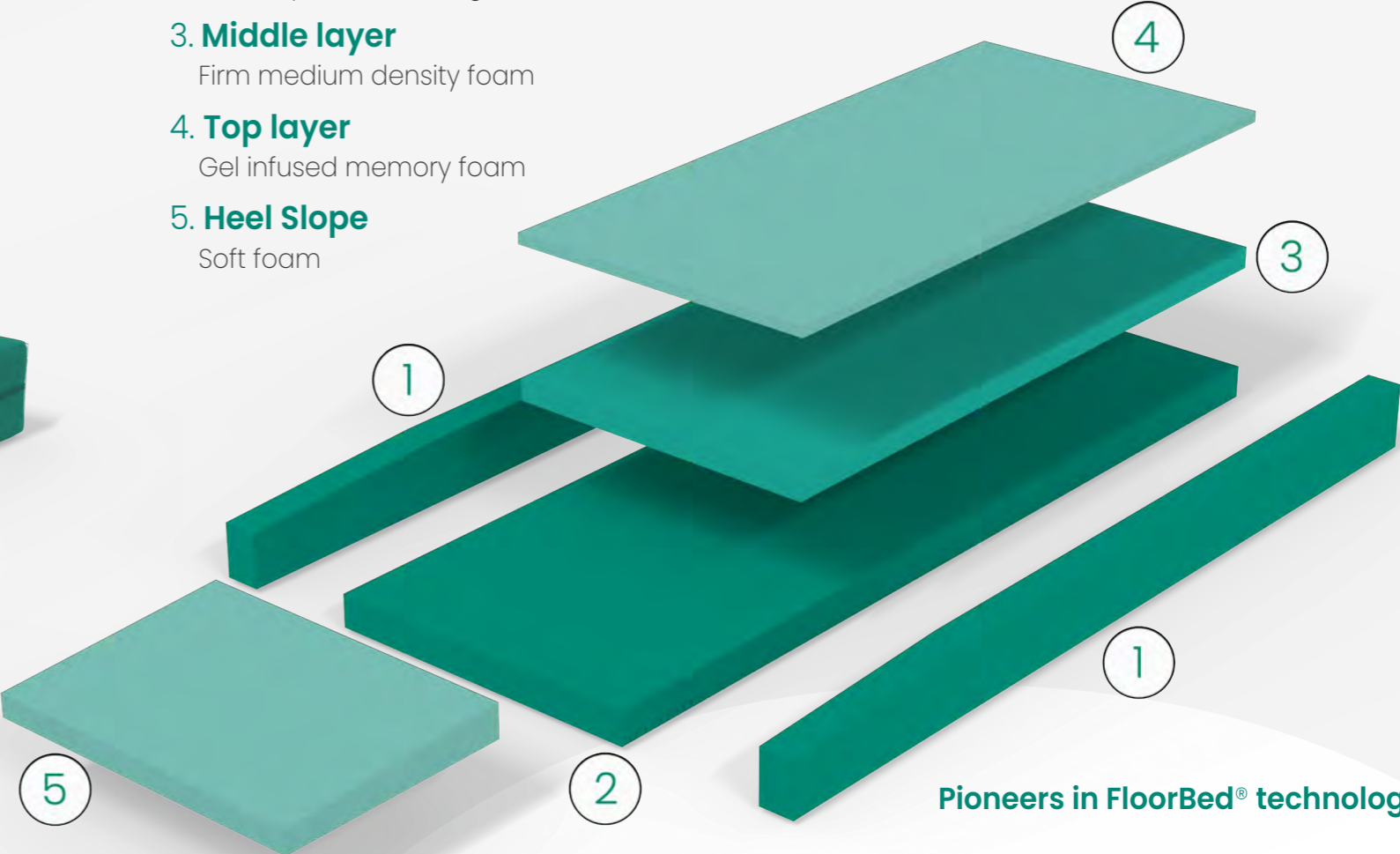

Pioneers in FloorBed® technology