



Mesh Safety Enclosure

The BBG canopy is a unique solution to a long standing concern by caregivers who keep asking for taller and taller safety sides. There are several engineering issues once you exceed a certain height including things like the unit being top heavy, no room to open the doors, etc.

For patients who are agile enough to need sides taller than our standard hi-side bed we recommend they consider adding a BBG mesh safety enclosure.



Key Points


- The unit is mesh with nylon reinforcement stretched over a tubular aluminum frame.
- Metal frame is bolted to corner posts.
- Adds approximately 15 inches interior height
- Mesh is white, fabric border comes in several colors
- Nylon sleeves cover metal frame
- Fabric covering is removable and washable.
- It is compatible with all BBG models



2 Mesh Choices

Standard

Heavy Duty (HD)

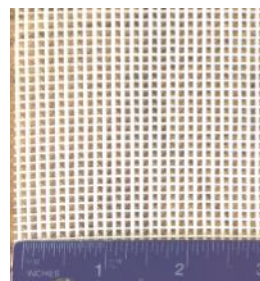

Order Info:

This canopy needs be ordered with the bed. Currently it CANNOT be retrofit.

This unit will use the Straight style head and footboards

This unit requires a modified squared off corner post.

This canopy is not compatible with any other brand bed

Suggested HCPCS: E0316

- ◇ Greater tear strength
- ◇ More visibility and airflow

Note: larger holes (1/4" dia) may not be suitable for tiny fingers

Fabric Colors



Polyester - Urethane Coated
Washable, mildew and water repellent

All images on this page are for general reference only. Electronic reproductions should not be used for exact color match. In addition, slight color variations will occur with any dyed fabric. We reserve the right to make product changes without notice.

availability of color
is not guaranteed



White


Updated January
2019



Grape




Beige

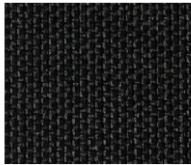


Black




Royal Blue

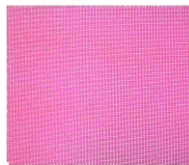

Hot Pink

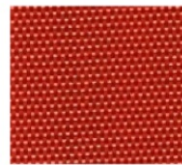

Red

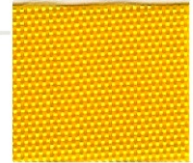

Yellow




Backwoods Camo




Navy Blue



Forest Green

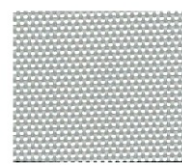

Silver Gray

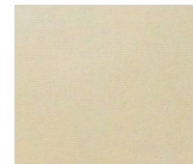

Ivory