

## Silver Sport Full-Reclining Wheelchair Single Axle



1P<sub>ACK</sub>


300  
⚙️

350  
⚙️

HCPCS  
K0001 / E1226

Item #	Seat Width	Armrest	Product Front Riggings	Color	UOM
SSP16RBDDA	16"	Detachable Desk Arm	Swing-away, Elevating Leg rests	Silver Vein	1/cs
SSP16RBDFA	16"	Detachable Full Arm	Swing-away, Elevating Leg rests	Silver Vein	1/cs
SSP18RBDDA	18"	Detachable Desk Arm	Swing-away, Elevating Leg rests	Silver Vein	1/cs
SSP18RBDFA	18"	Detachable Full Arm	Swing-away, Elevating Leg rests	Silver Vein	1/cs
SSP20RBDDA	20"	Detachable Desk Arm	Swing-away, Elevating Leg rests	Silver Vein	1/cs
SSP20RBDFA	20"	Detachable Full Arm	Swing-away, Elevating Leg rests	Silver Vein	1/cs

- Hydraulic reclining mechanism allows for infinite adjustments up to 180°.
- Arms have support extension for reclining position.
- Nylon upholstery is comfortable, lightweight and easy to clean.
- Carbon Steel frame with silver vein finish is attractive, and easy to maintain. (Figure A)
- Composite Mag-style wheels are lightweight and maintenance free.
- 8" front casters.
- Padded armrests provide added patient comfort.
- Wheels set back on frame prevents tipping.
- Precision sealed wheel bearings in front and rear ensure long lasting performance and reliability.
- Rear anti-tippers standard.
- Headrest extension with cushioned head immobilizer standard.
- Comes standard with swing-away elevating leg rests.
- Metal foot plates.
- Comes with push to lock wheel locks.

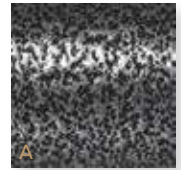

**specifications**

Width Open with 16" Seat	23"	Product Weight with 18" Seat	60 lbs.
Width Open with 18" Seat	25"	Product Weight with 20" Seat	61 lbs.
Width Open with 20" Seat	28"	Carton Shipping Weight with 16" Seat	66 lbs.
Width Closed	12"	Carton Shipping Weight with 18" Seat	67 lbs.
Height	50"	Carton Shipping Weight with 20" Seat	67.5 lbs.
Seat-to-Floor Height	20.5"	Weight Limit with 16" Seat	300 lbs.
Back Height	21"   31"	Weight Limit with 18" Seat	300 lbs.
Seat Depth	16"	Weight Limit with 20" Seat	350 lbs.
Length with Riggings	50"		
Product Weight with 16" Seat	59 lbs.		

**accessories**

<b>Anti Tippers</b>		<b>Anti Fold Bar</b>		
STDS4Y4712	Without Wheels	1/pr		
			STDS806 Anti-Fold Bar	
			1/ea	
<b>Seat Belts</b>		<b>Heel Loops</b>		
STDS850	Seat Belts, Auto-Clasp Type Closure	1/ea		
STDS851	Seat Belts, Velcro® Type Closure	1/ea		
STDS855	Seat Belts, Auto-Clasp Type Closure, Bariatric	1/ea		
STDS856	Seat Belts, Velcro® Type Closure, Bariatric	1/ea		
			STDS831 Heel Loops	
			1/pr	
			<b>Wheelchair Cane/Crutch Holder</b>	
			STDS1034	Wheelchair Cane/Crutch Holder
				1/ea