

Etac Swift shower stool/chair

Swift is a prize-winning series of shower stools/chairs designed with functionality and simplicity in mind. The basic shower stool can easily be transformed to a complete shower chair, by adding arm and back supports.

Stable on any floor

Swift has a unique design that adapts to uneven floors thanks to the flexible seat. The soft high quality ferrules adapt to the floor and ensures a good grip. All these features contribute to the users safety.

Changes with the circumstances

Swift is adaptable to different needs. It is possible to start using the shower stool and add arm supports and back support if the circumstances changes. Swift is easy to adjust to the required height with the telescopic legs. It can be set one step lower at the front for a slightly forward tilted seat to facilitate standing up. Soft pads which are specially suited for showering are available should extra comfort be needed.



See the movie about Swift on

Easy to handle

Swift is easy to handle thanks to its light weight. It comes flat-packed and is assembled without the need for any tools. Swift is also corrosion free.



Without arm supports, the attachments holes function as firm grips.

The arm supports have a stable fixation to ensure stability.

The back support has a practical grip for moving the chair



Easy to handle

The flat-packed Swift is easy to handle.



Safe surface

The structure of the surface ensures that Swift doesn't become slippery when wet.



Swift is also available in a low version.
Seat height 34.5–42 cm (13.6–16.5")



Standard seat height
42–57 cm (16.5–22.4")

Accessories



Back support
Green: 81701120
Grey: 84005069
Blue: 84005071



Back support pad

Green: 81701130
Grey: 84005072
The back support pad doesn't become slippery when wet and gives an instantly warm feeling. Material: Cellular polyethylene. Cleaning at max. 60°C (140°F)



Soap dish
Green: 81701140
Grey: 80209266
Blue :81701340



Seat pad

Green: 81701132
Grey: 84005073
The seat pad doesn't become slippery when wet and gives an instantly warm feeling. Material: Cellular polyethylene. Cleaning at max. 60°C (140°F)



Arm supports pair
Green: 81701110
Grey: 84005074
Blue: 84005075

Description

	Green	Grey	Blue
Swift shower stool	81701010	81701410	81701510
Swift shower chair with arm supports	81701030	81701430	81701530
Swift shower stool low	n/a	81701450	n/a
Swift shower chair low with arm supports	n/a	81701460	n/a

Dimensions

Seat height: 42–57 cm (16½–22½")
Seat height Swift low: 34.5–42 cm (13½–16½")
Total height with arm supports: 80–95 cm (31½–37½")
Total height Swift Low with arm supports: 72.5–80 cm (28½–31½")
Seat area: 54x41cm (21¼ – 16¼")
Total width, stool: 54 cm (21¼")
Total width, with arm supports: 54.5 cm (21")
Depth: 50 cm (19¾")
Depth Swift Low: 46.5 cm (18")
Width between arm supports: 45 cm (17¾")


Material

Polypropylene, aluminium

Weight

Shower stool: 3.1 kg (6.8 lbs)
Shower stool low: 2.9 kg (6.4 lbs)
Shower chair: 4.7 kg (10.4 lbs)
Shower chair low: 4.5 kg (8.8 lbs)

Colours

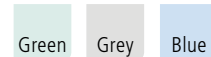

reddot



Clean the product with a solvent-free cleaning agent (pH level 6–9), or with a 70% disinfectant solution. Can be decontaminated at max. 85°C / 185°F if so required.