

The Levanto is a safe and comfortable standing and transfer device designed to help patients with some trunk balance and weakened leg muscles to stand up with a natural motion pattern. It features clip sling attachments and its compact design makes it the ideal option for smaller spaces.



Optimal Features

The Levanto is an electric sit-to-stand device that encourages a natural standing motion and is designed for standing, sitting, dressing, toileting and balance training for users with low core or leg strength. The shin pad angle adjusts with the patient while they are being raised for optimal comfort and security and the legs are easily spread electronically to allow direct access while seated (for example, in a wheelchair). The non-slip surface footplate offers stable support and is removable for easier cleaning and hygiene control.



Ideal for Home Care Use

The Levanto features a compact footprint that allows it to easily pass through narrow

doorways. It features a lifting arm with additional handles for comfort and easy usage. The shorter legs allow for easier rotation, especially in a home environment that may have less space for manoeuvring. With a LINAK battery and actuator, you are ensured years of trouble-free use.

Quality

Our products are designed and manufactured to meet all the requirements of relevant ISO, UL, IEC and EN standards. At every stage of the development process, every single product is given extensive attention in design, development and manufacturing. This means that a user of this product feels confident with its reliability and the caregiver feels confident with its safety.

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone: +46 8 665 35 00
Fax: +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street
Port Melbourne, VIC 3207
Phone: +61 1300 132 736
Fax: +61 (03) 9645 4954
info.au@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Accessories



Back Belts

Assortment of sizes and designs to accommodate most every lifting need.



Battery

Complete system with control box enclosed in a discreet module. Included with lift.



Hand Control

Slow or fast setting for arm up/down, electrical leg spreading. Included with lift.



Wheels

Rear wheels feature kick-brake. Diameter 125mm.

Lifting Situations



Use to/from a bed

Use to raise/lower a patient from/to a bed.



Use to/from a chair

Use to raise/lower a patient from/to a chair.



Use to/from the toilet

Use to raise/lower a patient from/to a toilet.

Features




Control Box

Emergency stop, emergency lowering and soft-start functions.



Lifting Arm

Lift arm extends 80cm / 31.5" to accommodate taller patients.

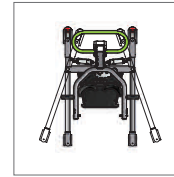

Spreadable Base

Electrical adjustment.






Tested to all relevant standards, including ISO 10535, hoists for the transfer of disabled persons.

Technical Specifications

Levanto			
Part No.	20000170	IP Rating	IPX4
Width	68.7cm / 27"	Turning Radius	128cm / 50.3"
Length	110cm / 43.3"	Front Wheels	10cm/4"
		Rear Wheels with Brake	12.5cm/5"
Height	Arm Lowered: 88.9cm / 35" Arm Raised ¹ : 164.8cm / 64.9"	Battery	24V DC/2.9Ah, 2x 12V/2.9Ah
Height from Floor	11.9cm / 4.7"	Weight	49kg / 108lb
Adjustable Leg [Min - Max]	Outside: 67.9cm - 113.1cm / 26.7" - 44.5"	Weight Capacity	160kg / 352.7lb
Sound Level	48dB (loaded) while lifting 52dB (loaded) while spreading legs		

¹ Measurement based on arm mast at highest setting

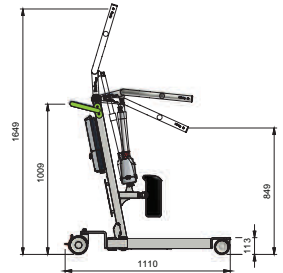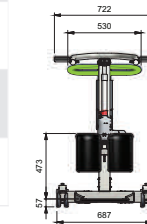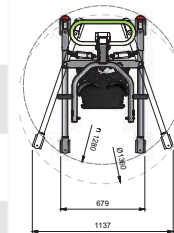

This product is manufactured in accordance with all relevant ISO standards.



Learn More About This Product



Follow Us on Social Media



facebook.com/humancarecda



instagram.com/humancaregroup



linkedin.com/company/human-care-ab



pinterest.com/humancaregroup



youtube.com/channel/UC7fyjNH3zTjRxOP1IDf2Da0