



INVACARE® MICROAIR® MA600  
THERAPEUTIC SUPPORT SURFACE



Yes, you can.®

**microAIR®**

# Invacare® microAIR® MA600

Alternating Pressure Low Air Loss Mattress System

Medium- High Risk  
Stage I-IV

3-in-1 Alternating Pressure System with  
Low Air Loss designed to help prevent and  
treat pressure injuries



**Highly vapor permeable and oversized pliable quilted nylon top cover** designed to provide low shear, friction and moisture protection



**Cable management system** keeps cords away from the ground to help prevent trip hazards




**Nylon non-slip base cover** designed to help reduce shear and friction during bed articulation




# Features & Benefits

**8 adjustable comfort settings** depending on the patient's therapy requirement

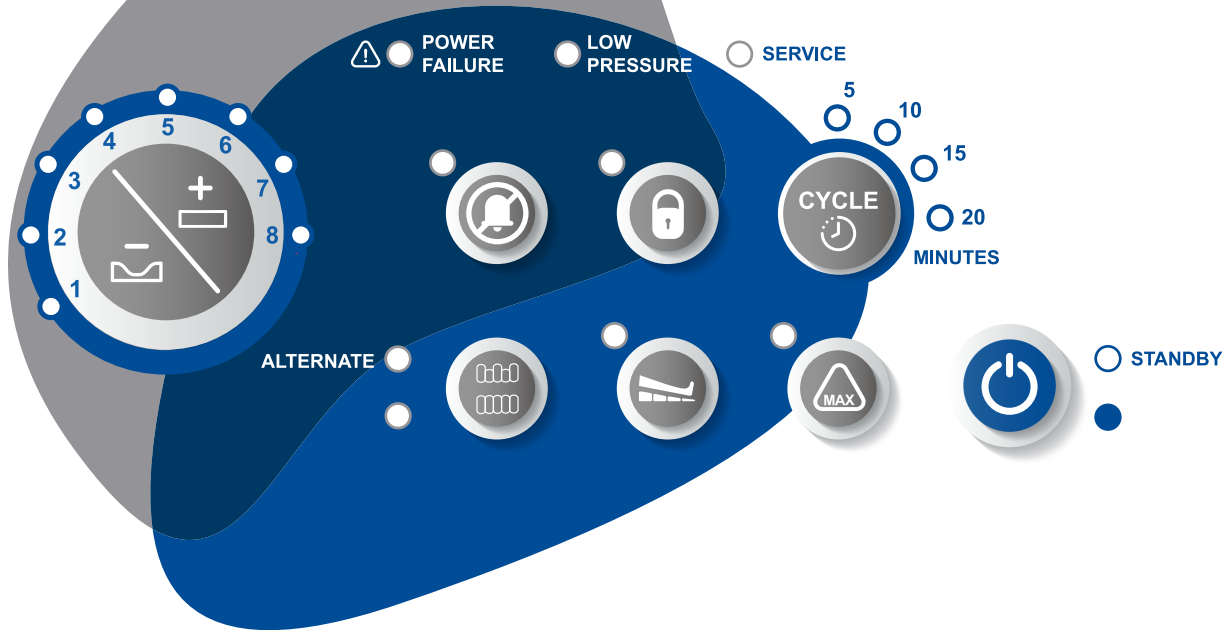

**Advanced 3-in-1 alternating function** inflates and deflates cells in a 3:1 cycle, providing active pressure relief by supporting 2/3 of the body at any given time

**Max Inflate Mode** provides quick maximum inflation within seconds to aid with transfers and nursing procedures

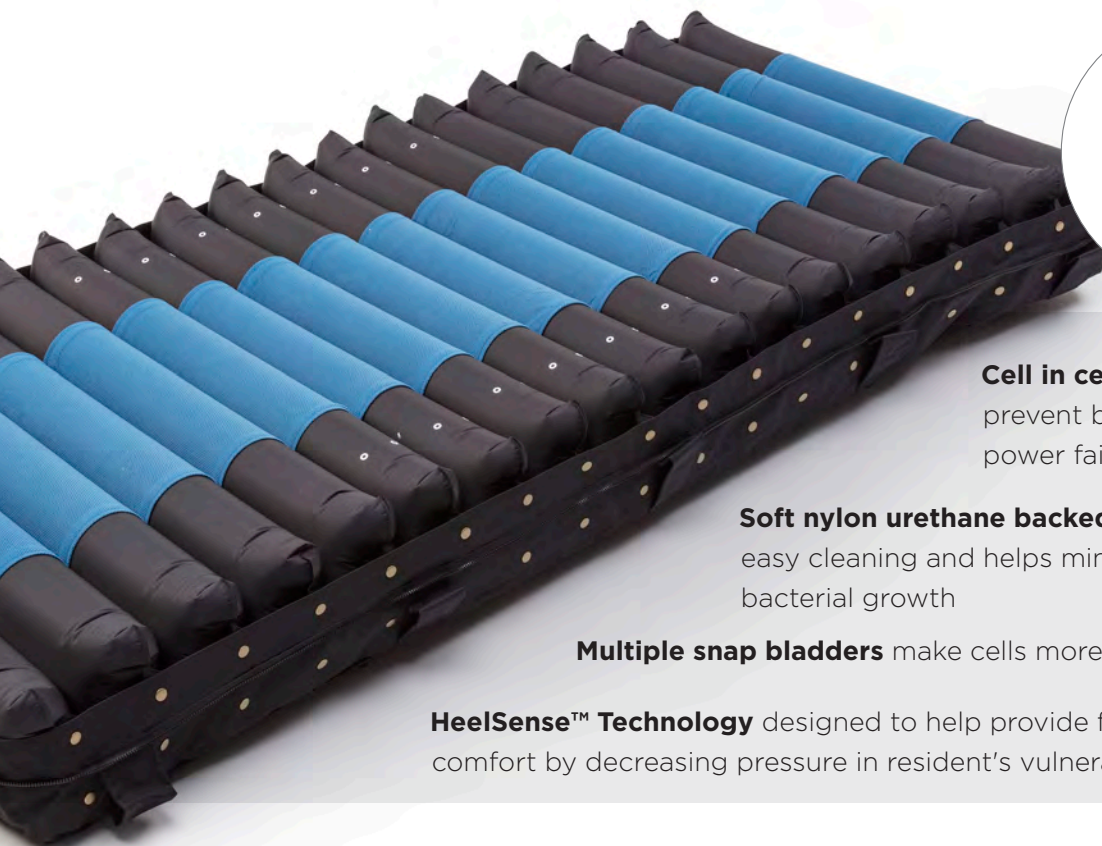

**Cell in cell design** is designed to help prevent bottoming out in the event of a power failure


**Soft nylon urethane backed cells** designed to provide easy cleaning and helps minimize the proliferation of bacterial growth

**Multiple snap bladders** make cells more secure with resident movement

**HeelSense™ Technology** designed to help provide further therapy and comfort by decreasing pressure in resident's vulnerable heel area

# Request a Quote or Place an Order Today!

## microAIR MA600 Alternating Pressure Low Air Loss Mattress System

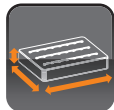

Model #	Dimensions	Weight Capacity	Flow Rate (Pump)	Therapy Type	Alarms	Limited Warranty
MA600	80" L 36" W 8" H	500 lb.	10 LPM	Alternating Pressure Low Air Loss	Audio & Visual	1 Year (Cover/Mattress) 1 Year (Power Unit)
MA600M	microAIR MA600 Mattress only					
MA600P	microAIR MA600 Pump only					

\*HCPCS Code: E0277

\*All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.