


INVACARE®  
**9000 XT**  
SERIES WHEELCHAIRS

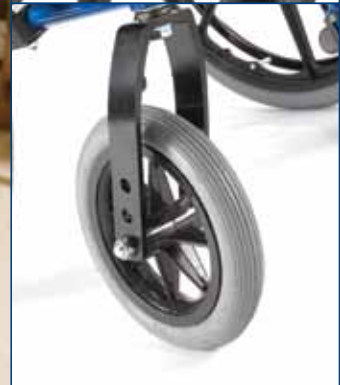

The **Invacare® 9000 XT Wheelchair** is a high performance, low maintenance and low total lifetime cost product. The 9000 XT wheelchair is simply the most versatile lightweight wheelchair available from Invacare, boasting an extensive offering of frame styles, seat-to-floor heights, back styles, widths, depths, frame colors and a wide range of options and accessories. The 9000 XT wheelchair components are interchangeable with all other 9000 series and Invacare® Tracer® series chairs.

The **Invacare® 9000 XT Recliner Wheelchair** is the number one choice for individuals who demand a truly versatile, yet moderately priced, lightweight wheelchair. The 9000 XT Recliner wheelchair can be configured into an endless number of seat widths, depths and heights, and has multiple frame and upholstery options. Components and accessories are interchangeable with the other 9000 series chairs.

### Invacare® 9000 XT Wheelchair Features

Dual axle positions allow variation of seat-to-floor heights  
 Offset axle position allows adjustment of center-of-gravity to increase maneuverability or stability to meet each user's skill level  
 Urethane rear tires mounted on "no flex" wheels offer superior performance  
 Urethane casters for better ride and lighter weight

### Invacare® 9000 XT Recliner Wheelchair Features

Dynamic recline range from 95° to 175°, +/- 5°  
 Rear axle base moves, providing stable center-of-gravity in recline position  
 Multiple wheel offerings allow variation of seat-to-floor heights  
 Urethane rear tires mounted on "no flex" wheels offer superior performance  
 Urethane casters for better ride and lighter weight  
 Low-maintenance, dual-sealed precision bearings  
 Heavy-duty inner liners keep seat and back from stretching  
 Stocked with attendant (standard) or patient-operated wheel locks

### HCPCS Code: K0004 & E1226 (recliner)



- 1. 10" Removable headrest
- 2. Dual axle positions
- 3. Urethane casters
- 4. Seating options available
- 5. Builder options available

Model #	Overall Width (2" increments)	Overall Height	Overall Depth w/o riggings	Folded Width	Seat Width (1" increments)	Seat Depth	Seat-to-Floor Height (1" increments)
9XT	21", 22", 23"-29"	36"	28"	12"	14"-20", 22"	16"-18"	13.5"-19.5"
9RC	23"-31", 28"	52"	32"	17"	14"-22" (2" inc)	16"-18"	17.5"-19.5"

	Back Height	Caster Height	Wheel Height	Weight Capacity	Shipping Weight	Carton Dimensions	Limited Warranty
9XT	15"-19"	6", 8"	20", 22", 24"	250 lb.	39 lb.	32"x13"x 37"	5 Years
9RC	34" w/headrest	8"	20", 22", 24"	250 lb.	42 lb.	32"x13"x 37"	5 Years

