

JAY[®]

BASIC FAMILY

- ≡ Introducing the newest addition to JAY's general use cushion offering – the high-performance, convenient JAY BasicPRO™



JAY® BASIC & BASIC^{PRO}™

≡ **The JAY Basic and Basic^{PRO} Cushions** are comfortable and durable foam cushions designed for clients at low risk of skin breakdown requiring mild stability. The Basic^{PRO} features a contoured base and easy to clean cover for performance and convenience. The Basic is comprised of a soft foam base and moisture-resistant cover for user comfort.

JAY Basic^{PRO}

Contoured Foam Base


- Provides positioning of the pelvis and thighs for increased comfort and support during use
- Durable, high-density foam provides superior performance throughout the life of the cushion

Easy to Clean

- Water-resistant Dartex® coated top cover surface is easy to wipe off and clean
- Optional Dartex®-coated Lycra® inner cover includes anti-wicking seam thread to protect the foam base from moisture build-up

JAY Basic

- Soft, mildly contoured foam provides the user a comfortable seated position for increased sitting tolerance
- Beveled base reduces “hammocking” effect associated with sling upholstery
- Moisture-resistant cover helps protect foam base

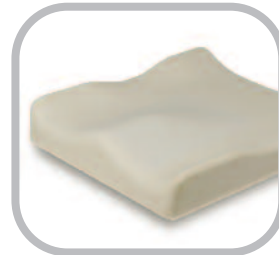

SPECIFICATIONS

JAY®

	Basic^{PRO}	Basic
HCPCS Code	E2601 & E2602	E2601 & E2602
Weight	1.9 lbs. (16" x 18")	1.9 lbs. (16" x 18")
Product Width	14" to 24"	14" to 24"
Product Depth	14" to 24"	16" to 20"
Product Height	3"	2"
Weight Capacity	14" to 21" widths: 300 lbs. 22" to 24" widths: 500 lbs	14" to 20" widths: 250 lbs. 22" to 24" widths: 300 lbs.
Outer Cover	Dartex®-coated Nylon Stretch and Ballistic Nylon	Moisture-Resistant