

JAY[®]

FUSION[™]

≡ Optimized for skin protection, stability and comfort



JAY FUSION™ CUSHION

≡ The JAY Fusion is an adjustable cushion designed for skin protection and features JAY Flow fluid or ROHO DRY FLOATATION™ air inserts in a package that's comfortable, stable and supportive. It incorporates the stability of the J3 cushion with the pressure management of J2 and it's easy to prescribe and order.

Simple to Fit

- Pelvic Loading Area (PLA) size based on anthropometric data and sized for cushion width
- PLA insert options (fluid or air) right sized for PLA

Simple to Order


- Smart part numbers
- Wide range of available modifications

Innovative Design Technology

- Soft upper base foam and steep well walls fit user comfortably and effectively
- Inner cover with Aquaguard moisture resistant zipper
- X-static® Anti-microbial outer cover with fluid option, or Stretch cover with air option

Complements JAY Adjustable Cushion Portfolio

- Stability from J3
- Pressure management from J2

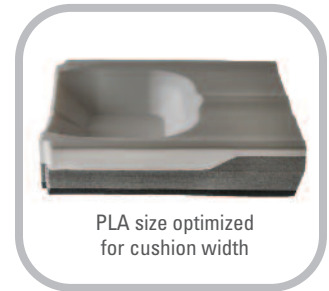

PLA size optimized for cushion width



Outer and inner covers



Anti-microbial silver thread

SPECIFICATIONS

JAY®

JAY Fusion Cushion	
HCPCS Code	E2622 & E2623
Weight	5.3 pounds (16" x 18" with fluid insert) 3.4 pounds (16" x 18" with air insert)
Product Width	14" – 24"
Product Depth	14" – 24"
Product Height	4"
Weight Capacity	14" – 21" widths: 300 lbs. 22" – 24" widths: 500 lbs.
Outer Cover	X-static™ (fluid), Stretch (air)