

Corner Seat & Nursery Table by Jenx

Corner Seat

An adjustable floor based seating system that provides support for young children from approx. 6 months to 10 years.

The Corner Seat includes an adjustable back height and can be used with or without the upper back support. An **Abduction Block** is also included to help keep the child positioned correctly; placing their hips neatly into the back of the seat.

Midline positioning is maintained by the corner design of the system and, if you prefer to have a flat back with wings at the side, there are small **Triangular Cushions** available to place at the back to provide this setting.

Nursery Table

A large and sturdy table that provides endless play opportunities.

The plastic activity surface can be used flat or angled as an easel and includes a height adjustable frame to meet the requirements of every child. It can also be used to great effect with the Corner Seat, with a short strap attaching the two items in use to prevent a child tipping the table over.

Surface area is easy to wipe clean and the cut out design to one side allows the child to get close to their table-top activities. Edging around the remainder of the surface also reduces the risk of objects sliding off the table.

Always Thoughtful **Outstanding features and benefits**



INNOVATIVE FEATURE
Trunk Positioning
The sturdy 'wings' and corner design help to maintain symmetry of a child's spine.

INNOVATIVE FEATURE
Floor Sitting
A simple first seat to introduce children to floor sitting. Long-leg sitting helps to maintain the length of the muscles behind the knee.


INNOVATIVE FEATURE
Height Adjustable
The Nursery Table can be easily adjusted for height. Adjustment is continuous rather than incremental, so every child can have their table at precisely the correct height for them. This feature allows the Nursery Table to be used by children using a Corner Seat or other floor seats, or by children who are able to sit unaided.

INNOVATIVE FEATURE
Shoulder Support
Adjustable upper support can be used to support the head or shoulders.

INNOVATIVE FEATURE
Sturdy and Bright
Sturdy table with a large, clear surface area; big enough for a child to play or work at most activities with ease.

INNOVATIVE FEATURE
Useful Cut Out Detail
The cut out on the table gives children the opportunity to be surrounded by the work surface, which means that even children who find it difficult to bring their hands forward will have their arms supported.

Corner Seat available in a range of 5 vibrant colours to choose from:

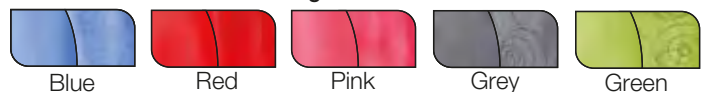

* Manufactured in a combination of both plain and patterned fabric

Always Configurable **Accessories**



Corner Seat Abduction Block

For abduction when long-leg sitting, the **Abduction Block** is adjusted to ensure the correct position for the child.



Corner Seat Triangular Cushions

For those who prefer to have a flat back with wings. The **Triangular Cushions** can also be used beside the legs to prevent external rotation of the legs.

Corner Seat Technical Specifications:

	Measurement	Units	Size 1	Size 2	Size 3
Age Range <small>(approx)</small>	approx.	mths - yrs	6 mths - 2 yrs	1 yrs - 6 yrs	4 yrs - 10 yrs
User Weight	max	kg	25	35	50
Total Back Support Height	min - max	mm	315 - 425	440 - 620	500 - 680
Lower Back Support Height	max	mm	165	235	235
Hip Width	max	mm	330	390	495
Seat Depth	max	mm	350	370	465

Side Lying Board Technical Specifications:

	Measurement	Units	Size 1	Size 2	Size 3
Table Size	length x width	mm	510 x 605	510 x 605	510 x 705
Table Height	min - max	mm	190 - 280	245 x 345	290 x 440
Table Angle <small>(From Horizontal)</small>	recline - forward	degrees	40° - 40°	40° - 40°	40° - 40°
Load Weight	max	kg	20	20	20
Distance Between Table Legs	width	mm	570	570	670

