

LECKEY®

Distributed by  Mobility

 Squiggles+

Stander





Standing Therapy, Reimagined.

The new Squiggles+ Stander offers additional features to the extremely versatile three-in-one Squiggles Stander for children ages 1-5 years. Now with hip abduction and improved function creating a superior clinical approach to standing therapy.

Posture



Hip Abduction

Hip Abduction now achievable, 60 degrees bilaterally. With a measurement gauge to monitor and aid repeatability. Hip laterals also shift in sync with the legs as they move through the arc of abduction.



Advanced Knee Support Design

Larger, more comfortable medial knee supports to distribute pressure, especially in abducted standing. Advanced wrap support to maintain knee in neutral flexion or extension.



Independent Footplates

Individually adjustable footplates allow you to set each footplate to different heights. This accommodates leg length discrepancies while providing support to both legs.

3-in-1

The Squiggles+ Stander is an extremely versatile 3-in-1 stander, offering prone, upright and supine in one product.



Prone

Provides anterior support and stability while encouraging weight bearing and use of upper extremities.



Supine

Support the back of the body and offer lateral and anterior supports to maintain alignment. Allows for partial weight bearing and promotes control of head, trunk, pelvis, knees and feet.



Upright

Designed for users with good balance and trunk control but cannot stand independently.

Function

Improved Pivot Base

Improved pivot base accommodates hip abduction and facilitates horizontal loading from 90 to 160 degrees. The new foot pedal eases caregiver effort while providing a more stable platform to transfer the child onto.




Activity Tray

Activity tray with improved height and adjustability accommodates the smallest users. Whether prone or supine, the tray with optional grab rail can be used to facilitate table top activities while in the standing position.



Height-Adjustable Push Handle

New height-adjustable push handle helps make positioning the child during standing therapy, easier. The height-adjustable push handle can be folded for ergonomic use and lowered for prone use.

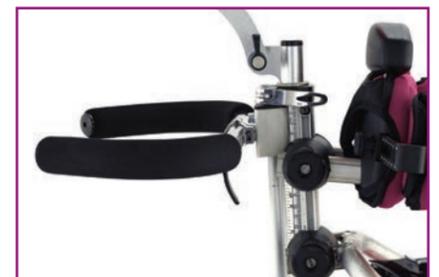

Squiggles+ Stander Specs

Age (approx)	1-5
Max User Weight	48.4 lbs.
User Height	29.5 in. - 43.5 in.
Chest Support Angle Adj.	+/- 30°
Chest Support Depth Adj.	2 in.
Chest Support Height (Footplate to Armpit)	21.5 in. - 32 in.
Chest Width (Distance between Laterals)	6.3 in. - 9 in.
Hip Support (Height from Footplate to Mid-Buttocks)	11.8 in. - 20.8 in.
Distance between Midline of Knees	5.5 in. - 8.4 in.
Footplate Angle	Plantarflexion 10° Dorsiflexion 10°

Support Frame & Pivot Base Specs

	Stander Support Frame	Base Footprint (Unfolded)	Base Footprint (Folded)
Weight	16.5 lbs.	14.3 lbs.	14.3 lbs.
Width	15.4 in.	23.2 in.	23.2 in.
Length	26.8 in.	37.4 in.	26.0 in.
Height	15.4 in.	20.3 in.	20.3 in.
Base Angle Range	90° - 160°	90° - 160°	90° - 160°

Tray Specs

	Prone	Supine
Tray Height (Footplate to top of Tray)	21.9 in. - 31.5 in.	18.5 in. - 31.5 in.
Tray Angle Adjustment	30°	40°

LECKEY[®]