

Canopy Enclosed Bed



34" wide or 42" wide

The Pedicraft Canopy Bed is now available in two widths.



27" High or 37" High

The Crib Height mattress level is now achieved by adding extensions to the legs. This raises the overall height by 10 inches. The extensions can be removed to convert the bed to a standard height.

#36850

All Pedicraft mattresses meet new Federal Safety Standard 16 CFR Part 1633 for the flammability of mattress sets.



Canopy™ Enclosed Bed

Offers the safety of an enclosed crib bed and a completely padded interior—no exposed metal parts. The Canopy enclosure is constructed of strong polyester mesh and Sure-Chek healthcare vinyl. The vinyl is bacteria resistant and the cover can be removed for cleaning. The Canopy Bed mattress is constructed of high density, fire-retardant foam covered in Sure-Check Comfort vinyl.

Canopy Bed includes:

- Padded Steel Bed Frame
- Mesh Enclosure
- 2-Layer Mattress
- Padded Bed Rails
- 6 Fitted Sheets
- Casters (2 locking, 2 non-locking)

The mesh panels unzip on all four sides for complete access to the patient. Locking zippers can be positioned where needed for the use of IV lines and feeding tubes while the bed is closed.

The Height Adjustable model provides convenience for the caregiver when a patient needs frequent medical attention. Raise the bed surface 10 inches by a hydraulic hand crank.



For extra protection the Pedicraft Canopy™ Enclosed Bed comes with padded bed rails, which are easily raised and lowered.


Options and Accessories



Padded Head and Foot Boards

For added support a head or foot board can be added.

#36848

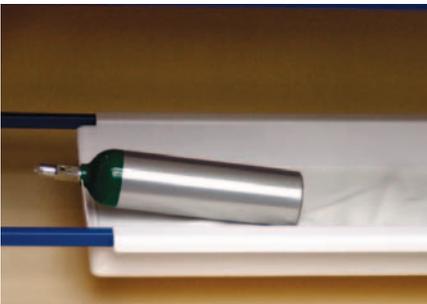

Storage Bin

Removable storage basket provides the perfect place to keep diapers, blankets and personal items. **#35822**



Head, Knee or Foot Elevation Wedge

Makes elevating the head, knees or feet quick and easy.

#36852

Optional Head Elevation Crank

Allows one end of the mattress to be elevated up to 30°, enabling the patient to sit up in bed.

Decorative Bed Skirts at No Extra Charge!

Give us an idea of your preferences and we'll choose the print to fit.

- A. PLANETS**
- B. HEARTS**
- C. SPORTS**
- D. PINWHEEL**

Sample swatches available.

Custom made skirts from fabric you purchase can be arranged.

