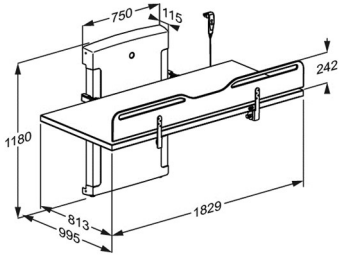




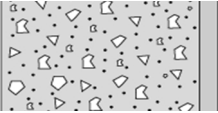

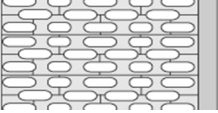
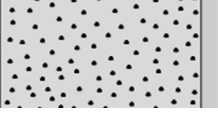
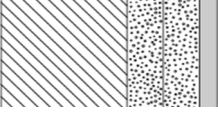
Product	Changing table, 32" x 72", manually foldable	
Item number	R8594572000	
	EAN:5708590358154	
	FDA:INQ/D124319	
Color	White	
Product description	<p>Changing table, 32" x 72", manually foldable. Electrically height adjustable with buttons placed in table front. Incl. safety rail, body belt and wired remote control. Excl. mounting kit. See Technical specifications (pdf) for information on mounting kits suitable for various wall types.</p>	
Adjustability	Height adjustable	
Adjustment height	700 mm	
Product capacity	<p>Maximum load: 440 lbs. When the table is loaded with 440 lbs, each mounting bolt supports 73 lbs. <i>This product has been tested statically with a load of 1100 lbs as the IEC 60601-standard requires testing with a safety factor.</i></p>	
Speed, upwards	<p>Speed up: 15 mm/s Speed down: 20 mm/s</p>	
Noise level	48 dB	
Mains connection (V)	110	
Power consumption (maximum load)	<p>Power consumption by maximum load: 70 W Power consumption by standby: 2 W</p>	
Declaration of conformity	<p>This product complies with the requirements of Council Directive 93/42/EEC concerning medical devices, as last amended by Council Directive 2007/47/EC of 5 September 2007. It also complies with all relevant requirements in the Council Directive 2006/42/EC (of 17 May 2006) on machinery and the Council Directive 2011/65/EC on RoHS.</p>	





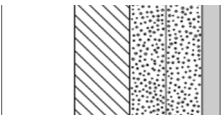
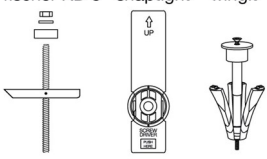
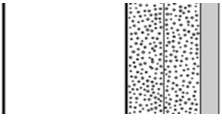
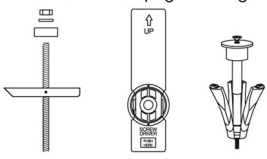
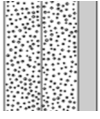

Product	Changing table, 32" x 72", manually foldable
----------------	--

Mounting in different wall types

	Concrete >B25 V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.
	Solid Brick >= Mz12 V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.
	Solid Sand Lime Brick >= Hlz12 V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.
	Aerated Concrete P4 V8658 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs GB14. To be ordered separately.
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. NB! Plugs are not needed. To be ordered separately.





Product	Changing table, 32" x 72", manually foldable	
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</p> <p>Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.</p>	<p>fischer KD 8 snaptight wingit</p> 
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch)</p> <p>Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.</p>	<p>fischer KD 8 snaptight wingit</p> 
	<p>Drywall (2 x 12,5 mm) water-resistant</p> <p>Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.</p>	

It is the responsibility of the installer to ensure that suitable wall fasteners are used.

