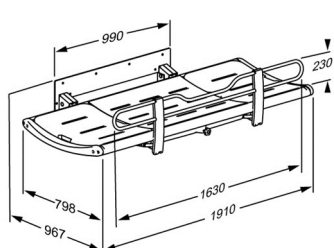




|  |  |   |
|--|--|---|
| <p><b>Product</b></p>  | <p>Shower and changing table 3000, fixed height</p>  |   |
| <p><b>Item number</b></p> <p><b>Color</b></p> <p><b>Product description</b></p>            | <p>R8435318021<br/> VVS:777543194<br/> EAN:5708590357478<br/> FDA:INT/D124320</p> <p>Lead grey</p> <p>Shower and changing table 3000, fixed height.<br/> Foldable. Incl. safety rail and water collection tray with integrated outlet.<br/> Surface length 70.9", outer measurements 39.6" x 75.2".<br/> Materials: Powder-lacquered stainless steel, aluminium, polyurethane foam and polystyrene.<br/> Excl. mounting kit.<br/> See Technical specifications (pdf) for information on mounting kits suitable for various wall types.<br/> <i>Replaces R8508 &amp; R8528.</i></p>   |  |
| <p><b>Adjustability</b></p> <p><b>Adjustment height</b></p> <p><b>Product capacity</b></p> | <p>Fixed mounting</p> <p>Maximum load: 440 lbs.<br/> When the table is loaded with 440 lbs, each mounting bolt supports 132 lbs.</p>   |   |
| <p><b>Materials</b></p> <p><b>Ergonomy</b></p> <p><b>Features</b></p>                      | <p>Aluminium, polyurethane foam, plastics, powder-coated steel.</p> <p>Pressalit offers a range of high quality, fixed or power adjusted changing tables designed to provide a stable platform suitable for changing adults and children in a variety of settings: schools, day care centers, institutions, hospitals and private homes.</p> <p>The showering and changing table 3000 series provides maximum safety and comfort. Besides an award winning design, the series includes comfortable and durable surfaces, including polyurethane foam and aluminium.</p> <p>A water collection tray with integrated water outlet prevents water from splashing on caregiver or the floor. The flexible hose attached to the center mounted outlet of the water collection tray, can be connected to either an outlet in the floor or on the wall.</p> <p>Adjustable neck support which can be moved to either end of the 70.9" changing tables.</p> <p>The slightly curved form makes the table extremely comfortable and safe.</p> <p>The safety rail can be folded easily up or down.</p> <p>The integrated pneumatic cylinders enable the table to be folded out of the way when not in use to maximize floor space.</p> |   |





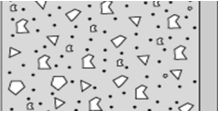

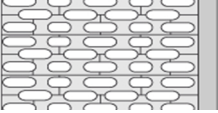
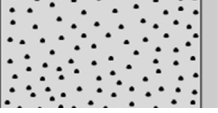
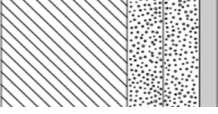
|                                  |   |
|----------------------------------|---|
| <b>Product</b>                   | Shower and changing table 3000, fixed height  |
| <b>Warranty</b>                  | Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is installed in an atmosphere heavy in salt or chloride, the warranty of the nursing bench will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product. The damaging effects can, however, be reduced if the product is showered with tap water on a daily basis.  |
| <b>Cleaning</b>                  | <p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b><br/>CaviCide</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b><br/>Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p> |
| <b>Declaration of conformity</b> | This product complies with the requirements of Council Directive 93/42/EEC concerning medical devices, as last amended by Council Directive 2007/47/EC of 5 September 2007.   |





|                |  |
|----------------|--|
| <b>Product</b> | Shower and changing table 3000, fixed height |
|----------------|--|

## Mounting in different wall types

|   |   |
|---|---|
|    | <b>Concrete &gt;B25</b><br>V8613 - Mounting kit, 8 lag bolts Ø10 x 70, plugs SX12, with grey covers.<br>To be ordered separately.   |
|    | <b>Solid Brick &gt;= Mz12</b><br>V8613 - Mounting kit, 8 lag bolts Ø10 x 70, plugs SX12, with grey covers.<br>To be ordered separately.   |
|  | <b>Solid Sand Lime Brick &gt;= Hlz12</b><br>V8613 - Mounting kit, 8 lag bolts Ø10 x 70, plugs SX12, with grey covers.<br>To be ordered separately.  |
|  | <b>Aerated Concrete P4</b><br>V8663 - Mounting kit, 8 lag bolts Ø10 x 70, plugs GB14, with grey covers.<br>To be ordered separately.  |
|  | <b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</b><br>V8613 - Mounting kit, 8 lag bolts Ø10 x 70, plugs SX12, with grey covers. NB! Plugs are not needed.<br>To be ordered separately. |





|                       |  |  |
|-----------------------|--|--|
| <p><b>Product</b></p> | <p>Shower and changing table 3000, fixed height</p>  |  |
|                       | <p><b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</b></p> <p>Wall fasteners (10 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.</p>           | <p>fischer KD 8    snaptight    wingit</p> |
|                       | <p><b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch)</b></p> <p>Wall fasteners (10 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.</p> | <p>fischer KD 8    snaptight    wingit</p> |
|                       | <p><b>Drywall (2 x 12,5 mm) water-resistant</b></p> <p>Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.</p>  |  |

It is the responsibility of the installer to ensure that suitable wall fasteners are used.

