

# UP-RITE I, II, III

## *Mid-Line Positioning Standing Frame*

The UP-RITE STANDER is a therapist designed upright standing frame which allows unobstructed access to the client from all sides supplying the ability for proper midline positioning for proper weight bearing and activities involving the upper extremities.



Shown with optional wheeled caster base.

### **THE UP-RITE STANDER INCLUDES:**

**All Cotton Body Support Straps** made from long-wearing, blue denim, breathable material.

**Tool Free Adjustments** for quick and easy midline fitting and support.

**Dual Side Posts** for complete access to client from anterior and posterior positions.

**Lexan Tray** with infinite height and depth settings encouraging upper extremity function, play and fine motor activities.

**Padded U-Shaped Knee Supports** available in three standard sizes include easy height, width and depth adjustment. *Optional Size-Adjustable Knee Pads offer additional width adjustment.*

**Adjustable, Padded Heel Support** with an open, “free movement” design allows for ease of foot placement. *Optional Shoe holders are available for additional positioning when needed.*

**Marine Grade Plywood Base** easily adapts for custom foot modifications. *Optional wheeled caster base is available to make the UP-RITE mobile.*


**Durable, All Steel Square Tube** construction for extra strength.

**Three Sizes** cover a wide range of clients.



# The Art of Standing

**UP-RITE I**



**Specifications**

Overall stander and max client height: 36"

Base dimensions:  
20" wide x 20" deep

Tray dimensions:  
18" wide x 14 1/2" deep

Knee pad cutout:  
3" wide x 2 3/4" deep

Unit weight: 21lbs.

**UP-RITE II**



**Specifications**

Overall stander and max client height: 48"

Base dimensions:  
24" wide x 24" deep

Tray dimensions:  
21" wide x 17 1/2" deep

Knee pad cutout:  
3" wide x 2 3/4" deep

Unit weight: 27 lbs.

**UP-RITE III**



**Specifications**

Overall stander and max client height: 60"

Base dimensions:  
32" wide x 32" deep

Tray dimensions:  
24" wide x 21" deep

Knee pad cutout:  
3 1/2" wide x 2 3/4" deep

Unit weight: 44 lbs.

UP-RITE standard features include: Height and depth adjustable clear table top; adjustable heel block; depth, height and width adjustable knee pads; one double and one single body support strap.

**O  
p  
t  
i  
o  
n  
s**

Chest and Hip Lateral Supports



Shoe Holders




Size Adjustable Knee Pads



V-Shaped Knee pads



Wheeled Caster Base

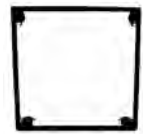

**Other Options:**

Double Body Support Strap  
Knee Straps

Foot Straps  
Head Strap



**PRIME ENGINEERING**  
A Division of Axiom Industries Inc.