



Mobile Shower Commode Chairs

Raz ... nothing else
comes close!



The Raz Advantage

Seats

The seat is the most important component of a Mobile Shower Commode Chair (MSCC). Raz has the largest selection of seats offered by any MSCC manufacturer: we offer 22 off-the-shelf models as well as custom.



IPAS (Ischial Pelvic Alignment System)

Raz seats are easily removable for cleaning and interchangeability. Our proprietary mounting brackets allow our seats to be adjusted fore / aft within a 2" range. These adjustments are critical for optimal alignment of the ITs within the aperture (commode opening). (see [Accessories Brochure](#))

Frame Sizes and Configurations

• AT: 5 widths, 2 lengths and 2 height ranges	• AT600: 5 widths
• AP/SP: 4 widths, 2 lengths and 2 height ranges	• AP/SP600: 5 widths
• CAT: 3 widths	
• ART: 1 width, 1 length and 2 height ranges	
• Züm: 1 width, 1 length	



Super-Low AT (5" lower)



Options and Accessories

Raz has the largest selection of options and accessories to create the optimal MSCC configuration. (see [Accessories Brochure](#))

Modular Construction

Our design philosophy centers around the principle of modularity. Raz MSCCs can be configured to accommodate the client's medical condition and functional needs.



SP Axle Plate



34"W AP600



Sacral Relief

Custom Solutions

When standard configurations don't meet specific needs, Raz will engineer and fabricate unique solutions including seats, frames and accessories.

Raz-AP Attendant Propel



- 4 widths and 2 depths
- Unrestricted front and side access
- 4 back heights (15" / 17" / 20" / 24")
- Adjustable back angle (85° / 90° / 95°)
- Tool-less height adjustment range (20"-24")
- 4 dual-locking casters
- Uses same frame as Raz-SP
- Stainless steel frame
- Optional reclining frame



Recessed, dipped front cross tube allows for better toilet bowl access and facilitates standing-pivot transfers



Unrestricted Side Access




Weight Capacity

Raz-SP Self Propel



Adjustable Axle Plate

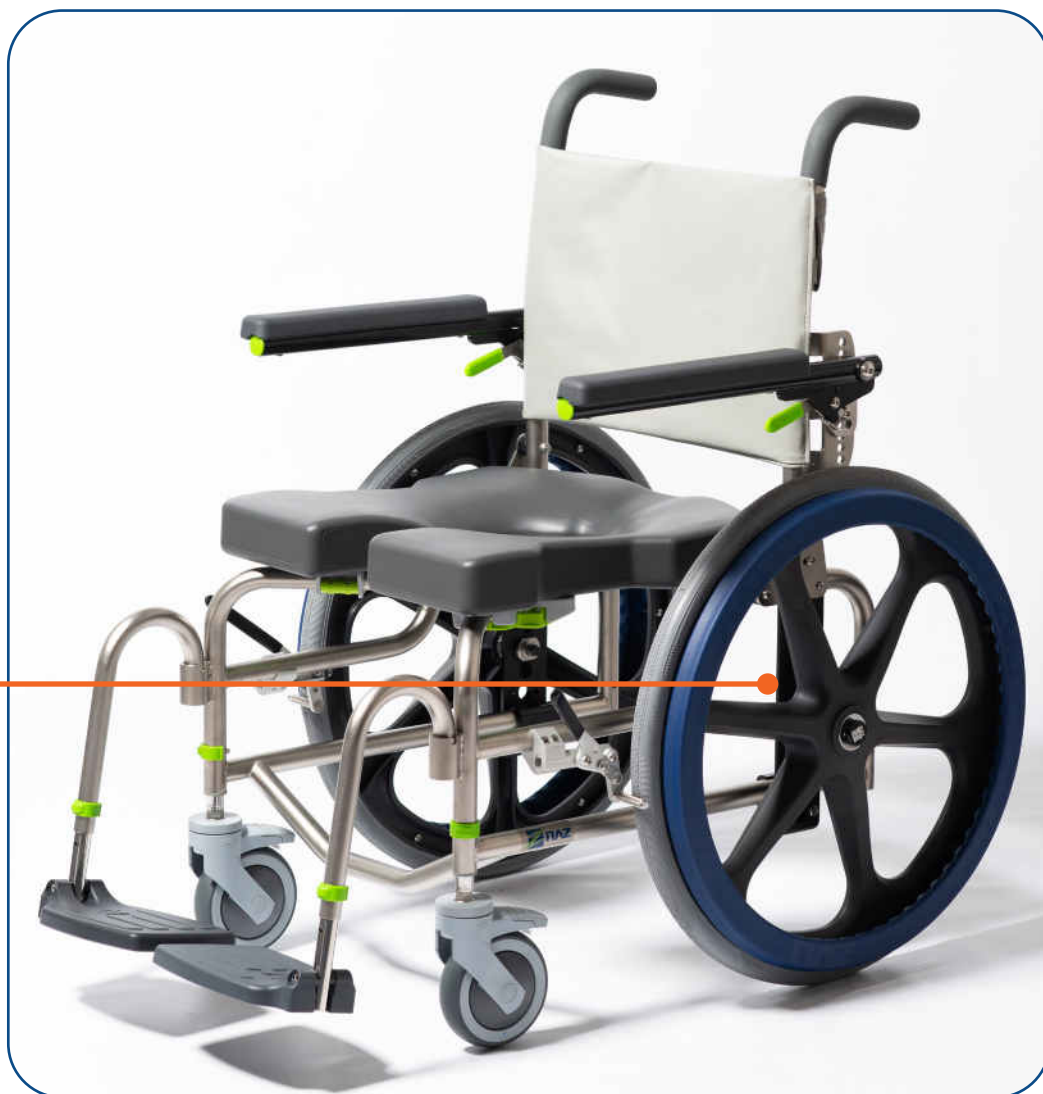

Ergonomic Handrim
with thumb groove and finger notches

- Center-of-mass adjustable axle plate (3" range)
- 22" or 24" wheels with ergonomic handrims
- 4 widths and 2 depths
- Unrestricted front and side access
- 4 back heights (15" / 17" / 20" / 24")
- Adjustable back angle (85° / 90° / 95°)
- Height adjustment range (20"-24")
- Super-Low frame available (16"-18")
- Uses same frame as Raz-AP
- Stainless steel frame
- Reclining Back optional

