



# Mobile Shower Commode Chairs

Raz ... nothing else  
comes close!



# The Raz Advantage

## Seats

The seat is the most important component of a Mobile Shower Commode Chair (MSCC). Raz has the largest selection of seats offered by any MSCC manufacturer: we offer 22 off-the-shelf models as well as custom.



### IPAS (Ischial Pelvic Alignment System)

Raz seats are easily removable for cleaning and interchangeability. Our proprietary mounting brackets allow our seats to be adjusted fore / aft within a 2" range. These adjustments are critical for optimal alignment of the ITs within the aperture (commode opening). (see [Accessories Brochure](#))

## Frame Sizes and Configurations

• <b>AT:</b>	5 widths, 2 lengths and 2 height ranges	• <b>AT600:</b>	5 widths
• <b>AP/SP:</b>	4 widths, 2 lengths and 2 height ranges	• <b>AP/SP600:</b>	5 widths
• <b>CAT:</b>	3 widths		
• <b>ART:</b>	1 width, 1 length and 2 height ranges		
• <b>Züm:</b>	1 width, 1 length		



Super-Low AT (5" lower)



### Options and Accessories

Raz has the largest selection of options and accessories to create the optimal MSCC configuration. (see [Accessories Brochure](#))

## Modular Construction

Our design philosophy centers around the principle of modularity. Raz MSCCs can be configured to accommodate the client's medical condition and functional needs.



SP Axle Plate



34"W AP600



Sacral Relief

### Custom Solutions

When standard configurations don't meet specific needs, Raz will engineer and fabricate unique solutions including seats, frames and accessories.

# Raz-AT Attendant Tilt



- 40° tilt range
- Gas struts help lift client out of tilt
- 5 widths and 2 depths
- 4 back heights (15" / 17" / 20" / 24")
- Adjustable back angle (85° / 90° / 95°)
- Height-adjustable frame (21"-25")
- Super-Low frame available (16"-20")
- 8° or 15° of anterior tilt
- Tilt Limiters available (no charge)
- Stainless steel frame



Pressurized Gas Struts

**400**  
lbs  
Weight Capacity

## Recline Module


- Available for Raz-AT / CAT / AP / SP
- 30° recline range (tilt limited to 25°)
- 2 reclining back heights (20" / 24")
- 2-Point Reclining Arms are standard
- Gas struts help lift client out of recline



# Raz-CAT Compact Attendant Tilt

## ... for growing kids and small adults!

- 3" shorter seat and base frames vs. Raz-AT
- Growable via modular back frames:
  - 3 back widths: 14" / 16" / 18"
  - 2 back heights: 15" / 18"
  - Seat depth adjusts from 11"-16"
- 40° tilt range
- Adjustable back angle (90° / 97°)
- Height-adjustable frame (21"-25")
- Seat with unique, tapered, oblong aperture
- Optional Recline Module
- Stainless steel frame

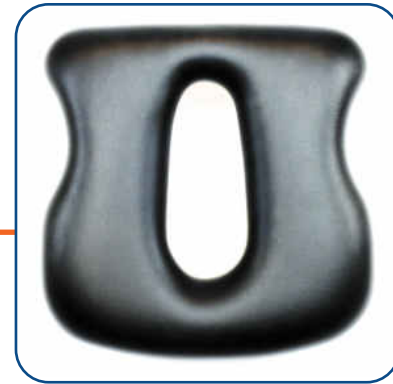

Standard CAT Seat

Tapered, oblong aperture (3.5"W at front and 4.5"W at rear). Aperture width increases as the seat depth increases.



3" shorter base frame length vs. Raz-AT



### Depth-Adjustable "Reach-Over" Back Frame

- Exclusive to Raz
- Allows for back frame to be narrower than seat frame
- Arm supports can be positioned optimally for trunk width of client
- Seat depth range 11"-16"
- 10" arm support lengths available for shorter seat depths (13" is standard)