



## *The ENHANCER<sup>®</sup> DRY FLOATATION<sup>®</sup> Cushion* by ROHO, Inc.

The ENHANCER DRY FLOATATION Cushion offers you a unique option in wheelchair seating. By offering a contoured form with a two manifold system, the client/user can benefit from positive midline channeling of the legs along with increased lateral stability. In addition, the special design of the ENHANCER cushion inhibits sacral sitting while enhancing tissue viability.

The ENHANCER's unique contoured shape is the result of flexible air cells of varying heights. The two manifold system allows you to: 1) adjust the interior bucket seat section for reduced pressure on the ischia and coccyx while 2) adjusting the two outer perimeter rows for enhanced stability and positioning. The two design features combined offer you a cushion that shapes to the user's body providing superior pressure distribution and minimal shear forces that are characteristic of every ROHO<sup>®</sup> DRY FLOATATION cushion, while offering you the added benefit of progressive position.




HIGH PROFILE<sup>®</sup>  
QUADRO<sup>®</sup>



nexus<sup>®</sup>

**CROWN<sup>®</sup>**  
THERAPEUTICS, INC.  
The Body's Partner

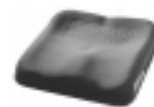

AirLITE<sup>™</sup>



LOW PROFILE<sup>®</sup>  
QUADRO<sup>®</sup>

# Features and Benefits

**Variable Cell Heights:** Provides for channeling of the legs and enhances lateral stability and body positioning.

**Washable, Two-Way Stretch Cover:** Allows for proper immersion into the cushion and “breathes” with every shift in position.

**Adjustable:** Shapes to the individual and to their specific seating needs.



**Ventilation:** Individual air cells allow for air circulation between the cells.

**Lightweight:** The ENHANCER is approximately four (4) pounds making it easy to transport and handle and minimizes fatigue of the individual.

**Durable Neoprene Rubber:** Flame resistant and resistant to urine, oils and medicines. Easily cleaned with soap or water and can be steam autoclaved. Product can be repaired extending the overall life of the cushion.

**Dual-Manifold Design:** Permits specific adjustments to allow for reduced pressures on bony prominences while independently promoting better positioning and increases stability.

**24 Month Warranty**

## Specifications

The ENHANCER comes complete with hand inflation pump, repair kit and operating instructions.

Construction—flame resistant neoprene rubber

Weight—approximately four (4) pounds

Cleaning—soap and water, or can be steam autoclaved

MODEL	CHAIR SIZE
ENH88C	fits 15"-16" widths & 15"-16" depths
ENH89C	fits 15"-16" widths & 16"-17" depths
ENH98C	fits 16"-17" widths & 15"-16" depths
ENH99C	fits 16"-17" widths & 16"-17" depths
ENH910C	fits 16"-17" widths & 18"-19" depths
ENH109C	fits 18"-19" widths & 16"-17" depths
ENH1010C	fits 18"-19" widths & 18"-19" depths

Can be custom made to fit a variety of seating surfaces and sizes.

## Applications

**Most effective for users:**

- ◆ at moderate risk for skin breakdown
- ◆ needing some degree of abduction/adduction
- ◆ with poor sitting balance
- ◆ who sacral sit

**CROWN®**  
THERAPEUTICS, INC.