

# Mobility at home

SAVARIA CEILING LIFTS



Home care

 savaria®



**Jonathan Benoit**  
Research & Development

“ We knew we could make a better lift by questioning the details and focusing on simplicity and reliability

### Peace of mind, reducing risks

Whether you are a patient or a caregiver, a Savaria ceiling lift provides security for lifting and transferring in a home setting. The lift helps prevent the risk of injury to caregivers, including back strains that are common occurrences with manual lifting. Patients can feel peace of mind and security when moving or transferring using the Savaria lift. For some, a Savaria ceiling lift is the solution for leaving the hospital sooner and recovering in the comfort of their home.

### Better technology, better lift

Like today's quality power tools and even electric cars, Savaria lifts use lithium-ion batteries with an optimized transmission to bring a host of benefits including:

- **More cycles on one charge:** 50% more than the industry average
- **Less charging time:** typically 15 minutes a day, or 2 hours if fully depleted
- **Power:** faster lifting speeds for less time spent in the lift
- **Long lasting:** batteries expected to last 5 years, no worrying about dead batteries and related service calls





### Portable lift + Gantry

Paired together, this system provides simple lifting from wheelchair to bed. The Savaria gantry design features low profile feet for easier wheelchair maneuvering in small spaces and the gantry is fast and easy to assemble. Add a reacher arm for easier positioning of the lift.



### SureClip system

The sling clip opens out rather than pushing inward. This detail creates a more foolproof system to attach the sling, eliminating errors that cause the sling to detach from the lift during use.



### Slings

Slings are offered in universal, hygienic and walking styles in 3 sizes. All slings feature a variety of adjustments to create a customized and comfortable fit. Toss them in the washing machine for easy laundering – no stiffening materials to remove first.

### Easy to use

The Savaria design team studied every detail in order to make improvements over traditional designs. Rather than adding extra bells and whistles, the team focused on making the lift perform reliably and intuitively. From the first time you try it, you will know how to use it with its straightforward design.

### Powerfully light, simply strong

Benefiting from lithium-ion battery technology and unencumbered design, Savaria ceiling lifts are light – the portable lift is only 11 lbs and the fixed (installed) lift weighs only 17 lbs.

### Fixed lift + FastTrack

For permanent installations, choose the Savaria fixed track in a custom configuration to fit your home including room covering or multi-room designs. The motorized option provides even more ease to move around.




## SPECIFICATIONS

Product weight	Portable: 11 lb (5 kg) Fixed model: 17 lb (7.5 kg) Fixed motorized model: 19 lb (8.5 kg)
Maximum load	Portable: 286 lb (130 kg) or 440 lb (200 kg) Fixed: 600 lb (272 kg), 440 lb (200 kg), or 286 lb (130 kg)
Noise level	Maximum 54 dBA
Lifting speed	2.2"/sec (5.5 cm/sec) at 0 lb 2"/sec (5 cm/sec) at 286 lb (130 kg) 1.8"/sec (4.5cm/sec) at 440 lb (200 kg)
Battery type	Lithium-ion, 25.2V 2500 mAh
Lifts per charge	Approximately 80 cycles of 24" (61 cm) at 200 lb (91 kg)
Battery charging	Full capacity in approximately 2 hours
Product life expectancy	10 yrs with only 1 battery replacement
Warranty	2 years
Compliance	Meets all international certification and compliance standards for medical devices

## OPTIONS & ACCESSORIES

- Motorized fixed lift in 3 capacities
- Portable lift and freestanding gantry in 2 capacities
- FastTrack system available for simple or extensive room covering designs in modular straight, 45 degree and 90 degree curve sections as well as XY trolley and turntable system
- Hardware for wood, concrete and steel structure installations
- Reacher arm for portable lift
- Universal sling in breathable mesh fabric with and without head support. Hygienic sling (toileting/bathing) and walking sling (rehab use) in easy-clean solid fabric. All slings available in small, medium, or large with a wide spectrum of adjustments for optimal fit and comfort.

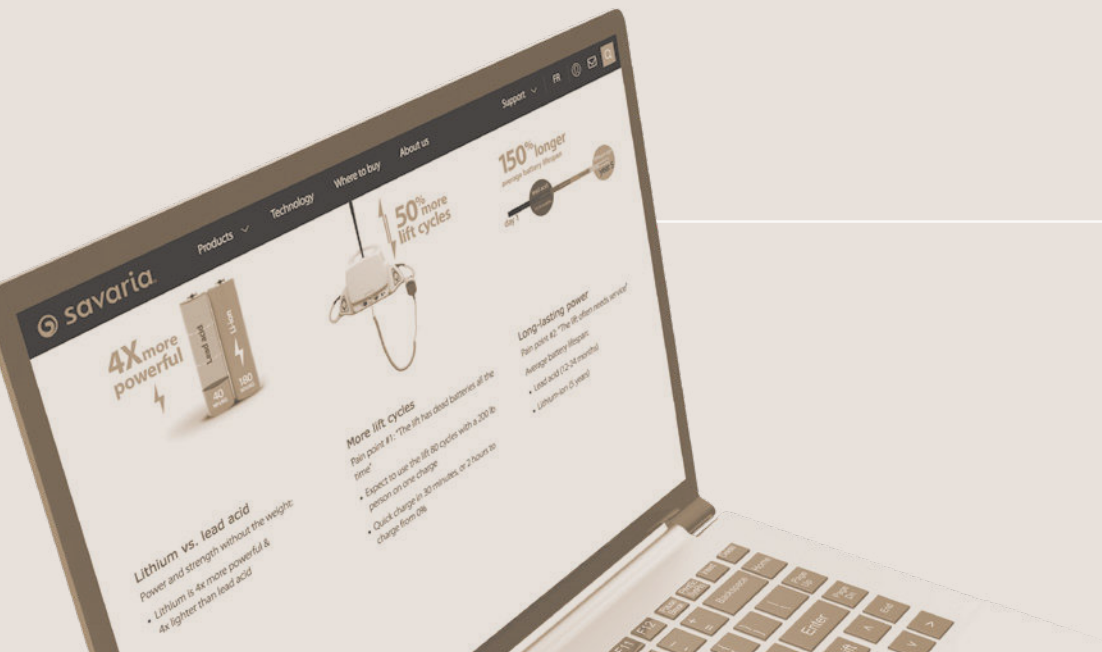

**ABOUT SAVARIA** Savaria Corporation is one of North America's leaders in the accessibility industry with wholly-owned manufacturing facilities in the US, Canada and China. Savaria offers a large portfolio of mobility products including home elevators, wheelchair lifts, stairlifts and accessible vans. The company also owns Span-America, makers of high quality PressureGuard® mattresses and medical beds. The company's mission is centered around mobility for home, public spaces and transportation, creating better mobility for life.