




Patent No. 7,434,282 B2  
AND OTHER PATENTS PENDING

First in the development of complete CELL LOCK technology, Star offers you the StarLock Cushion. This is our premier wheelchair cushion offering you the ultimate in aggressive positioning, unparalleled stability and the highest level of therapeutic prevention and treatment of pressure ulcers you've come to expect from a Star Cushion Product.




5" cell height

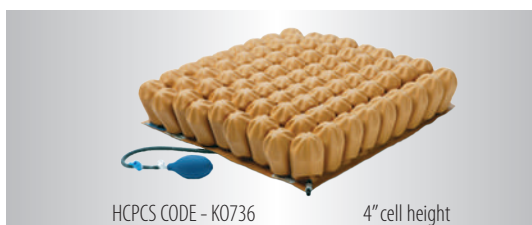


HCPCS CODE - K0736

4" cell height

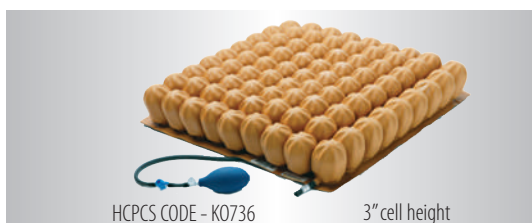

HCPCS CODE - K0736

3" cell height

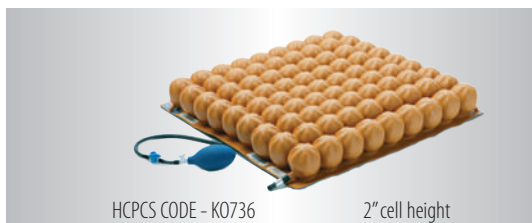

HCPCS CODE - K0736


2" cell height



4" cell height with cover

- Available in 5", 4", 3", and 2" cell heights to better accommodate your seating and positioning needs.
- Variable air flow control provides versatility when addressing hard to fit clients.
- LOCK TECHNOLOGY allows you to easily address asymmetrical seating issues.
- LOCK TECHNOLOGY allows you to isolate individual air cells and customize with no contact areas without the added expense of an actual custom cushion.
- The ability to isolate air flow between cells enhances stability and assists in transfers and overall posture or trunk control.
- Maintains skin integrity by providing a therapeutic, low pressure, low shear, low friction seating environment.
- Cost effective seating option for a wide range of seating and positioning needs.
- Light Weight, Approximately 3 pounds./1.36 kilograms.
- Available in Pediatric and Adult sizes.
- 24 Month Manufacturers Limited Warranty.
- Comes complete with fitted cover, hand pump and repair kit.
- Designed of high quality, flame resistant latex free neoprene rubber.
- No Weight Limit

## Skin Protection, Positioning and Stability

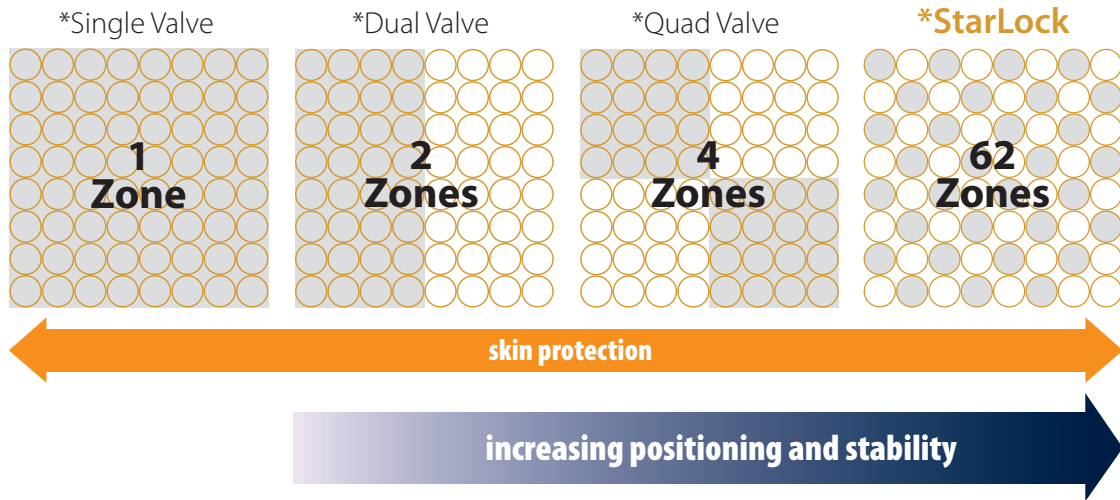

*\*Configuration based on a 64 celled air cushion. (8 x 8 array)*

Innovative Starlock Technology creates a multitude of zones offering unprecedented positioning and stability. The Starlock set-up is quite simple... adjust the Starlock like a traditional single valve air cushion. Using the attached bulb pump, inflate the locking mechanism to secure client position. When locked, the amount of air volume created in each individual air cell by virtue of the set-up procedure will remain...creating a very personalized seating surface, without the cost of a custom cushion.



**STARLOCK**<sup>®</sup>