

NEW



EASyS Modular S

More freedom of movement
for an active day


Thomashilfen

More head support, more comfort

This is **NEW** ...

“ The perfect choice for children who need more positioning support in combination with the best possible freedom of movement. The short back in addition to the multi-dimensional adjustable headrest ensures ideal adaptation and support for active children. A large variety of adjustments for optimal seating, positioning and comfort in motion are the basis for taking part in everyday activities. ”

Small height and depth adjustable arm supports are positioned close to the child's body, for relaxed sitting

Multi-dimensional adjustable headrest for optimum adaptation and support (four possibilities)



Choose between 4 different **HEADREST** possibilities: **Thomashilfen**



NECK SUPPORT:
supports the sensitive neck area during head movements



SHELL:
gives maximum support for the whole head



ADJUSTABLE:
both sides adjust for individual positioning



ADAPTER:
an option for a custom-made headrest



A single release lever on the headrest tubes is used for:

- seat tilt-in-space
- separating the seat from the base



Soft suspension absorbs vibrations and ensures a smooth ride without compromising steering

The different headrest possibilities and adjustments allow perfect head support for children while promoting an anatomically correct head and sitting position.

Short back and separate headrest provide freedom of movement in the torso and upper extremities. Children sit in a stable position while taking part in their environment.

Sophisticated ergonomics allow it to adapt as the child grows. All necessary therapeutic adjustments can be handled hassle-free without tools.



Easy seat assembly and rapid change of viewing direction – without tools:

- forward facing
- rear facing

140° resting position offers comfort and relief for the child's spine.



Much more storage space available under the seat:

- washable basket (inclusive)
- vent tray for medical devices (optional)

EASyS Modular S

More freedom of movement for an active day

Chassis & accessories (selection)



Q Chassis:
for kindergarten, school
and at home with
comfortable height and seat
tilt-in-space adjustment



Handy accessory plate
for medical devices, e.g.
oxygen tank or suction
equipment



Flexible lateral trunk
supports and a 4-point
pelvis belt for positioning
the upper body and the
pelvis




Grip rail and flexible
butterfly/chest harness for
a comfortable and secure
position of the upper body



Comfortable sleeping bag,
padded inside or with
woven fur for cold days



Abduction block for
guiding the thigh and
padded h-belt for holding
the upper body.

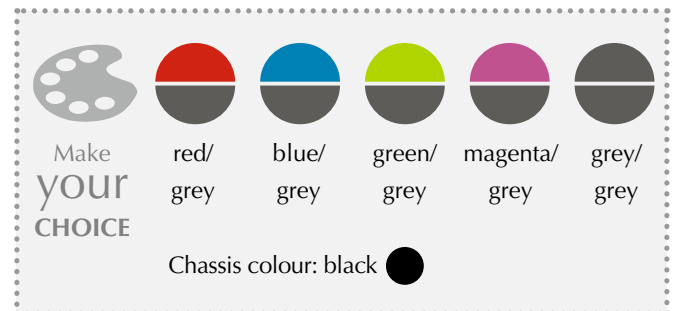

Choose between 5 cool
upholstery styles. The
cuddly cushions are
breathable and washable.

Technical data (seat)	EASyS Modular S 1	EASyS Modular S 2
Seat depth	19 - 30 cm 7.5 - 11.8"	27 - 38 cm 10.6 - 15"
Seat width	19 - 26 cm 7.5 - 10.2"	27 - 35 cm 10.6 - 13.8"
Seat height	47 cm / 18.5"	46 cm / 18.1"
Back height	28 - 42 cm 11.0 - 16.5"	34 - 53 cm 13.4 - 20.9"
Lower leg length	18,5 - 30,5 cm 7.3 - 12.0"	18,5 - 37 cm 7.3 - 14.6"
Hip angle	90° - 140°	90° - 140°
Knee angle	90° - 180°	90° - 180°
Foot angle	78° - 100°	78° - 100°
Turning range of seat	180°	180°
Size of hip guide (w x h)	23 x 12 cm 9.1 x 4.7"	29 x 11 cm 11.4 x 4.3"
Size of foot rest (w x d)	28 x 18 cm 11 x 7.1"	32 x 21 cm 12.6 - 8.3"
Weight seat unit	8,8 kg / 19.4 lbs	11 kg / 24.3 lbs
max. user weight seat unit	35 kg / 77.2 lbs	40 kg / 88.2 lbs

Technical data	A1 chassis	A2 chassis
Seat tilt-in-space driving direction	-10° to +35°	-10° to +35°
Seat tilt-in-space face-to-face	0° to +35°	0° to +35°
Push bar height	71.5 - 117 cm 28.2 - 46.1"	71.5 - 117 cm 28.2 - 46.1"
Wheel size swivel wheels (front)	18 cm / 7.1"	18 cm / 7.1"
Wheel size (rear)	25 cm / 10"	25 cm / 10"
Weight chassis	9.8 kg / 21.6 lbs	10.9 kg / 24 lbs
max. carrying capacity basket	5 kg / 11 lbs	5 kg / 11 lbs
max. carrying capacity (incl. seat)	50 kg / 110.2 lbs	55 kg / 121.3 lbs

Technical data stroller	EASyS Modular S 1	EASyS Modular S 2
Overall dimensions (l x w x h)	99,5 x 59,5 x 114 cm 39.2 x 23.4 x 44.9"	107 x 69 x 120 cm 42.1 x 27.2 x 47.2"
Folded dimensions (l x w x h)	85 x 59,5 x 60 cm 33.5 x 23.4 x 23.6"	81 x 69 x 69 cm 31.9 x 27.2 x 27.2"