

zippie®

IRIS

ROTATION-IN-SPACE
TECHNOLOGY DESIGNED
JUST FOR KIDS!

The Zippie IRIS® wheelchair offers a superior combination of positioning, maneuverability, and now increased portability with the new XLOCK® folding crossbrace option.



Check out the IRIS web page

FEATURES AND BENEFITS

Intelligent Rotation in Space

WITH PEDIATRIC-SPECIFIC DESIGN

Available with 40° or 55° of patented Intelligent Rotation in Space technology. Designed specifically for pediatric users to ensure a smooth tilting motion. Helps your child remain comfortable and secure while tilting and is easy for you to control.



Optional XLOCK® Folding Crossbrace

The patented folding crossbrace fully locks when unfolded for superior rigidity and folds into a compact size for convenient transport.



Built-in Growth


Offers up to 4" of built-in depth growth and 2" of built-in width growth with the XLOCK frame or optional width-adjustable strut tubes. Adjustable components, including armrests and footrests, help to maintain proper positioning during development.

WC-19 Transit Approved for Safe Bus Transit


WC-19 transit approved for use as an occupant's seat in a motor vehicle when used with an approved Wheelchair Tiedown and Occupant Restraint System in accordance with SAE (J2249).

Z-finity™ System for Lower Extremity Positioning

Lower leg length setting can be adjusted between -.5" and 21" from the seat and footrest angle can be adjusted infinitely to achieve the best positioning as your child grows.

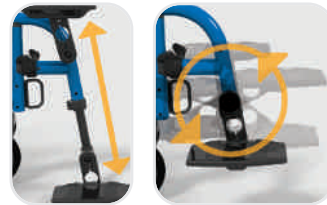

Options for Caregiver Convenience

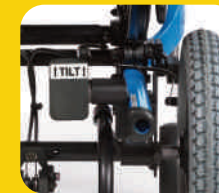

- Vertical mount vent tray provides convenient access to either LTV or Trilogy ventilators

- Foot Release Tilt Actuator allows for tilting without exposed cables (pictured left)

- Push-button hub locks located on the back canes are compact, easy to access, and easy to operate



- Angle Adjustable Stroller Handle with 260° of angle adjustability fits varying caregiver heights and enables pushing from a tilted position (pictured left)

- Rocker arm cover prevents debris from collecting on the hardware (available on rigid frame only)

HCPCS Code:

Folding: E1234
Rigid: E1233

Transit Option: Yes*

Tilt:

40°: -5 to 35°, 0 to 40°, +5 to 45°
55°: -5 to 50°, 0 to 55°, +5 to 60°

Chair Weight:

Folding: 35.8 lbs
Rigid: 32 lbs

Transport Weight: (w/o footrests & rear wheels)

Folding: 26.6 lbs
Rigid: 22.8 lbs

Weight Capacity: 225 lbs (w/o seating)

Overall Width: 17.5" to 25.5"

Seat Widths:

Folding: 10" to 16"
(grows to 18")
Rigid: 10" to 18"

Seat Depths: 10" to 18"

Front Seat-to-Floor Height:

13 to 19"
w/drop seat:
11" to 17.5"

Rear Seat-to-Floor Height:

13" to 19"
w/drop seat:
11" to 17.5"

True Fit or One-Time Free Growth Kit:

IRIS: Yes
IRIS SE: No

Transit Option: Yes

Conforms with ANSI/RESNA WC-19 standards. Refer to Quickie® & Zippie® WC-19 Transit Securement Supplement (MK-100072) for additional information. Supplement is available for download at www.sunrise-medical.com.

Aztec Gold	Black	Black Cherry	Black Opal
Blue Opal	Candy Blue	Candy Purple	Candy Red
Electric Blue	Evergreen	Glow	Green Apple
Hot Sparkle Pink	Mauve Pink	Rootbeer	Sparkle Silver
Sunrise Orange	Yellow	Titanium Color Paint	Matte Black
Matte Black Cherry	Matte Purple	Matte Electric Blue	Matte Evergreen
Mossy Oak Camo	Stars and Stripes	Zebra*	*Black Zebra Kolorfusion is available on any base color.