

zippie®

ZONE™

Where's your ZONE?



The **Zippie ZONE** is an ultra lightweight, option-rich, pediatric rigid chair that features flip-up footrests, sideguards with built-in hublock access and an abducted frame with inset backrest. The ZONE ergonomically fits some of the smallest children while still allowing the greatest wheel access for independent mobility.



Light and Efficient

Ultra lightweight:

- 10" X 10" configuration weighs as little as 14 lbs.

Compact Efficiency:

- 7000 series aluminum and Shapeloc Technology create strength while ensuring lowest possible weight

Growth and Adjustments

Growth Capability:

- 3" of built-in frame depth adjustment
- 2" seat sling depth adjustment
- Frame width growth is as simple as 4 tubes (seat, back, camber and footrest)

Simple Adjustability:

- Backrest and axle system adjust together and maintain center-of-gravity to backrest position
- Center-of-gravity adjusts independently in 1/4" increments up to +2.25"

Easy and Cool

User Accessibility:

- Flared seat and inset backrest to accommodate small kids and provide better access to the wheels for propulsion
- Wheellock in sideguard is easier to reach for rider, and one lever actuates both wheels
- Hublock on seat rail- release button positioned out of the way when pushing
- Flip-up angle adjustable footrest for easy transfers in and out of the chair

Cool Factor:

- 19 frame color options including Glow and Root Beer
- 7 anodized color package options
- Zippie spoke guard design




JAY's Innovative Pediatric Cushion!

The JAY® Zip™ was truly designed for skin protection and positioning through the use of anthropometrics of hip breadth, weight and ischial tuberosity spacing.

The Zip features a contoured foam base that's comfortable, effective and convenient.

Order with your Zone today!

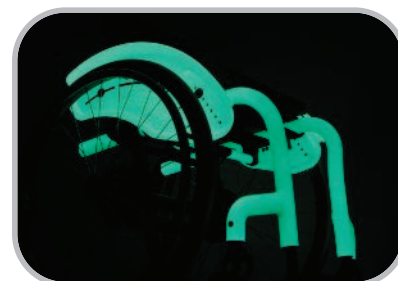

ZIPPIE ZONE SPECIFICATIONS

HCPCS Code:	E1235
Overall Weight 10" X 10":	14 lbs.
Weight Capacity:	165 lbs.
Seat Width:	8" - 16"
Seat Length:	8" - 18"
Center of Gravity:	- 0.25" - +2.25"
Camber:	0°, 3° or 6°
Back Height Adjustments:	6" - 8" up to 14" - 18"
Front Seat Height:	14"
Rear Seat Height:	13"
Frame Angle:	80°
Warranty:	Limited/Lifetime Frame